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An experience of more than 70 years at the manufacture of flasks and systems for storage and transport of liquid gases, tempering and insulation of laboratory devices made of glass or metal enables KGW - ISOTHERM to produce and distribute high quality products for science and industry. The periodically controlled quality of the materials used in manufacture and the quality control while production of the KGW - ISOTHERM products have been developed a high confidence between user and manufacturer. Modern methods in development (CAD) and manufacture of products such as a certified quality management have been increased the confidence in KGW - ISOTHERM products. Products other than the products described in the catalogue, KGW - ISOTHERM manufactures devices and containers according to the customers requirements in series or individual fabrication.

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Material - Glass

All glass devices produced by KGW-ISOTHERM are strictly made of Borosilicate-glass 3.3 DIN/ISO 3585.

Chemical Characteristics

Hydrolytic resistance:
according DIN-ISO 719 (98 C)
Acid resistance
according DIN-ISO 1776
Alkaline resistance class
according ISO 695-A2

Optical Characteristics

Spectral region, where the absorption is negligible low:
310 - 2200 nm

Physical Characteristics

linear expansion factor:
(at 20-300 C): 3,3 x 10-6 1/K
Density: 2,23 g/cm³
specific thermal capacity: 910J/kgK
Transformation temperature 525 C

Surrounding temperature Dewar vessels
+150°C to -200°C

Cleaning

The cleaning of KGW - ISOTHERM laboratory glass devices can be realised manual or in a laboratory dish washer. The commerce offers a large program of purifiers.

Important: Do not clean glass devices with scouring powders or scouring swabs such as other devices that could damage or hurt the glass surface.

Standards and Regulations

All glass devices are manufactured according the German "AD-Merkblätter", the EU Pressure Equipment Directive 97/23/EG, and the DIN EN ISO 16496 "Devices with vacuum insulation". If there are standards for joint components like spherical ground joints or conical ground joints, KGW-ISOTHERM take them to account. (e.g. DIN 12242-1 and DIN 12244-1)

Quality Control

Each KGW - ISOTHERM product is defeated by an individual quality control and designation such as documentation of all production data.

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Application

- flexible applications in laboratory
- cooling of small samples
- storage:
 - liquid nitrogen (LN2)
 - dry ice (CO₂)
 - other coolants
- as insulating flask

Advantages

- high insulating ability
- easy handling
- made of Borosilicate glass 3.3
- deliverable:
 - fully silvered (standard)
 - silvered with viewing stripes
 - shortened

Versions

- C: blue coated metal cover
G-C: blue coated metal cover with side grip

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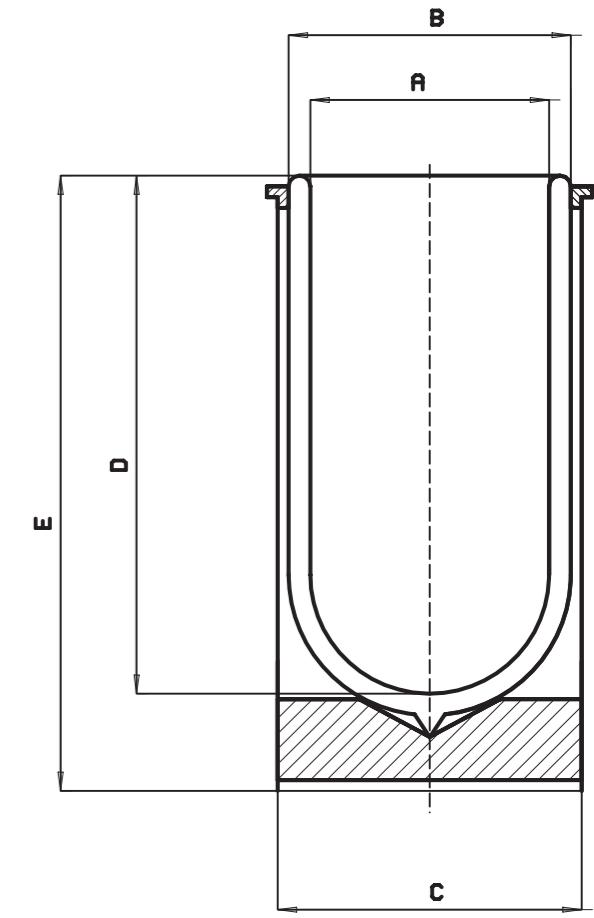
Application

- glass refill for Dewar flasks cylindrical, see page 6-7

Advantages

- high insulating ability
- easy handling
- made of Borosilicatglass 3.3
- deliverable:
- fully silvered (standard)
- not silvered
- silvered with viewing stripes
- shortened

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Typ F-C / Type F-C

Application

- in plants where a gas-tight joint is required
- suitable as insulating flask
- if using gassing coolants, take care of the max. over pressure of 0,1 bar !

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Versions



C: blue-coated metal cover
 A: glass refill (without cover)

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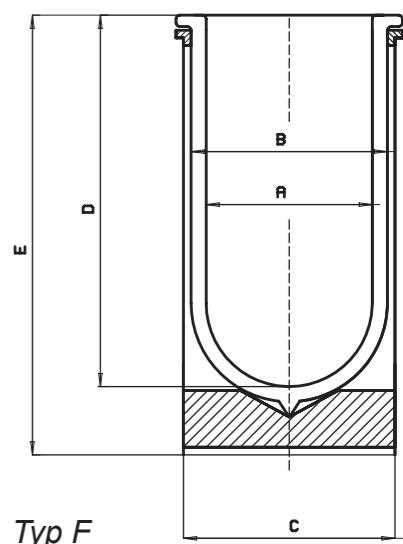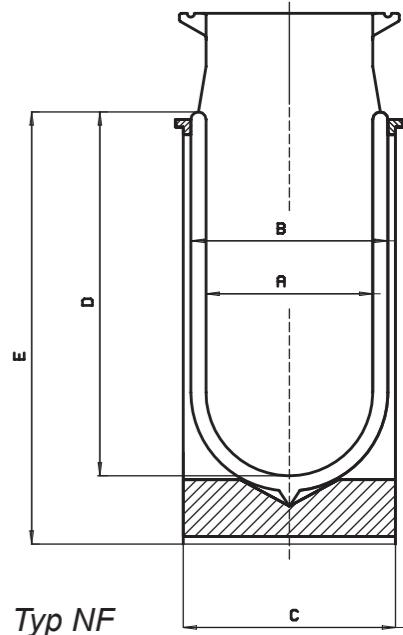Application

- suitable on magnetic stirrer
- flexible applications in laboratory
- cooling of small samples
- storage:
 - liquid nitrogen (LN2)
 - dry ice (CO_2)
 - other coolants
- suitable as insulating flask

Advantages

- low evaporation rate
- easy handling
- aluminium cover:
 - antimagnetic
 - stainless
- made of Borosilicateglass 3.3
- deliverable:
 - fully silvered (standard)
 - silvered with viewing stripes
 - shortened

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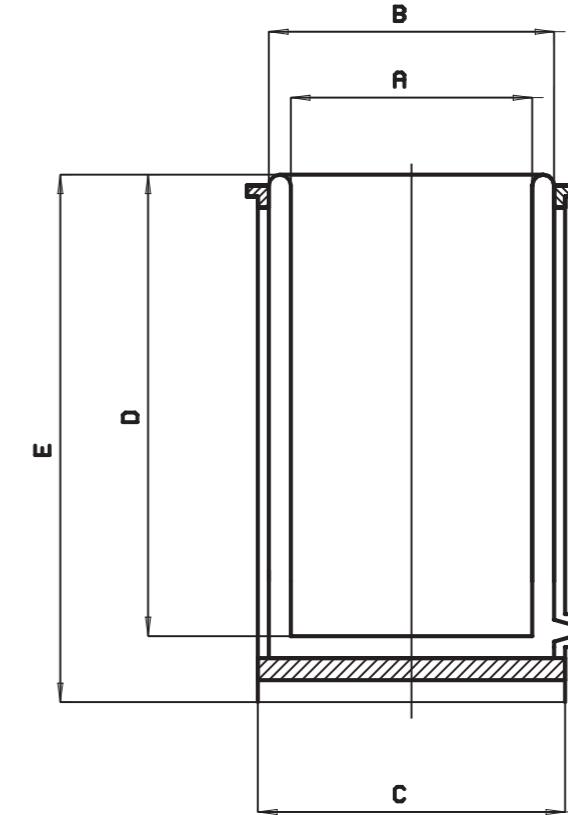
Typ SCH-A / Typ SCH-A

Application

- tempering of round bottom flasks
- cooling of small samples
- suitable as a hot or cold bath
- suitable on magnetic stirrer

Advantages

- easy handling
- aluminium cover:
 - antimagnetic
 - stainless
- made of Borosilicateglass 3.3
- fully silvered

Versions

CAL: structured aluminium cover
 A: glass refill (without cover)

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Typ B-E / Type B-E

Application

- storage and transport of temperature-sensitive goods e.g. dry ice (CO_2)
- flexible applications in laboratory
- cooling of small samples
- as high quality insulating flask
- when storing liquid nitrogen (LN2) a hole (\varnothing ca. 2 mm) in centre of the lid is necessary.

Advantages

- robust design
- made of Borosilicate glass 3.3
- with insulating lid held by two spring clips
- incl. handle made of aluminium

Versions

B: blue-coated metal cover
 B-E: stainless steel cover
 A: glass refill (without cover)

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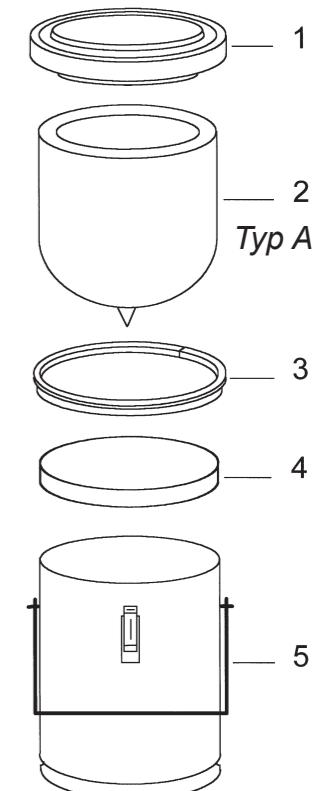
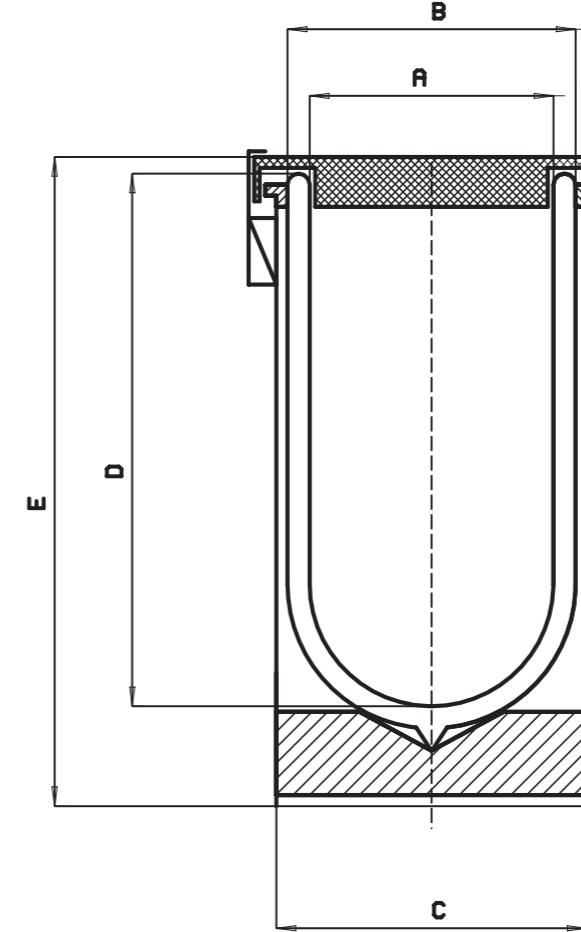
Typ / Type 27BE-Bio-V3

Area of application

- laboratory technique
- medical technique
- biotechnology
- storage and short internal transport of biological or pharmaceutical samples

Features of performance

- reliably and economically
- Dewar flasks according to DIN EN ISO 16496
- glass material according to ISO 3585
- nominal capacity from 2 to 4 litres
- protective casing out of stainless steel
- for storage in CO₂
- pressureless space for coolants

Description of complete flasks

Types and accessories

- B-E = protective casing out of stainless steel
- carrying handle = pure aluminium
- lid = polyethylene, white
- clamps = stainless steel, electrolytic polished
- V2A sample tubes for indirect cooling in CO₂

Temperature ranges

- glass refill from -200°C to +150°C,
- lid made of polyethylene
- Storage of samples at up to -196 °C

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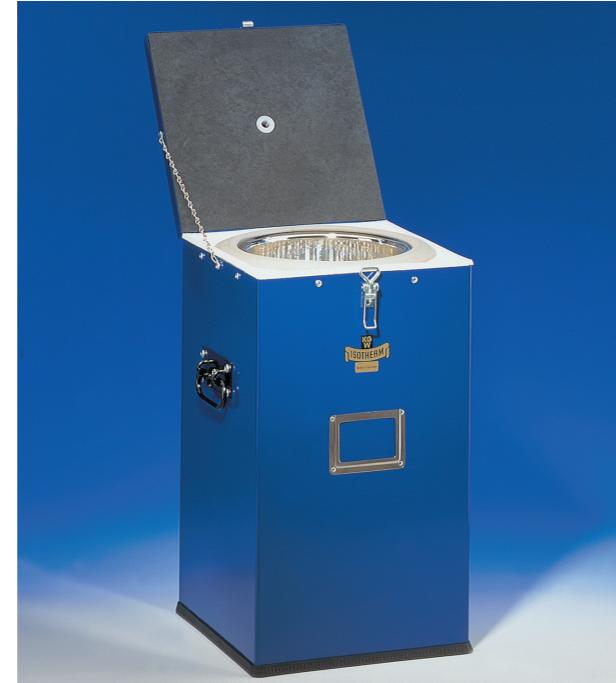
Typ 133 / Type 133

Application

- storage and transport of temperature sensitive goods e.g. dry ice (CO_2)
- suitable as insulating flask

Advantages

- good tipping stability
- robust design
- large aperture
- made of Borosilicateglass 3.3
- incl. - insulated lid
- quick release clamp
- two side grips

deliverable:

- fully silvered (standard)
- silvered with viewing stripes
- shortened

Versions

- blue-coated metal cover
- Type 43 to 46 are used as glass refills of type 131 to 134

Accessories

- linen bag for easy handling of dry ice (CO_2)
- insulating discs, recommendable for storing of dry ice
- metal roller base for type 133-134

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Application

- storage and transport of liquid nitrogen (LN2)

Advantages

- low evaporation rate
- good tipping stability
- made of Borosilicate glass 3.3
- incl. loose lid
- incl. handle
- deliverable:
 - fully silvered (standard)
 - silvered with viewing stripes

Versions

AL: aluminium cover, brushed
 A: glass refill (without cover)

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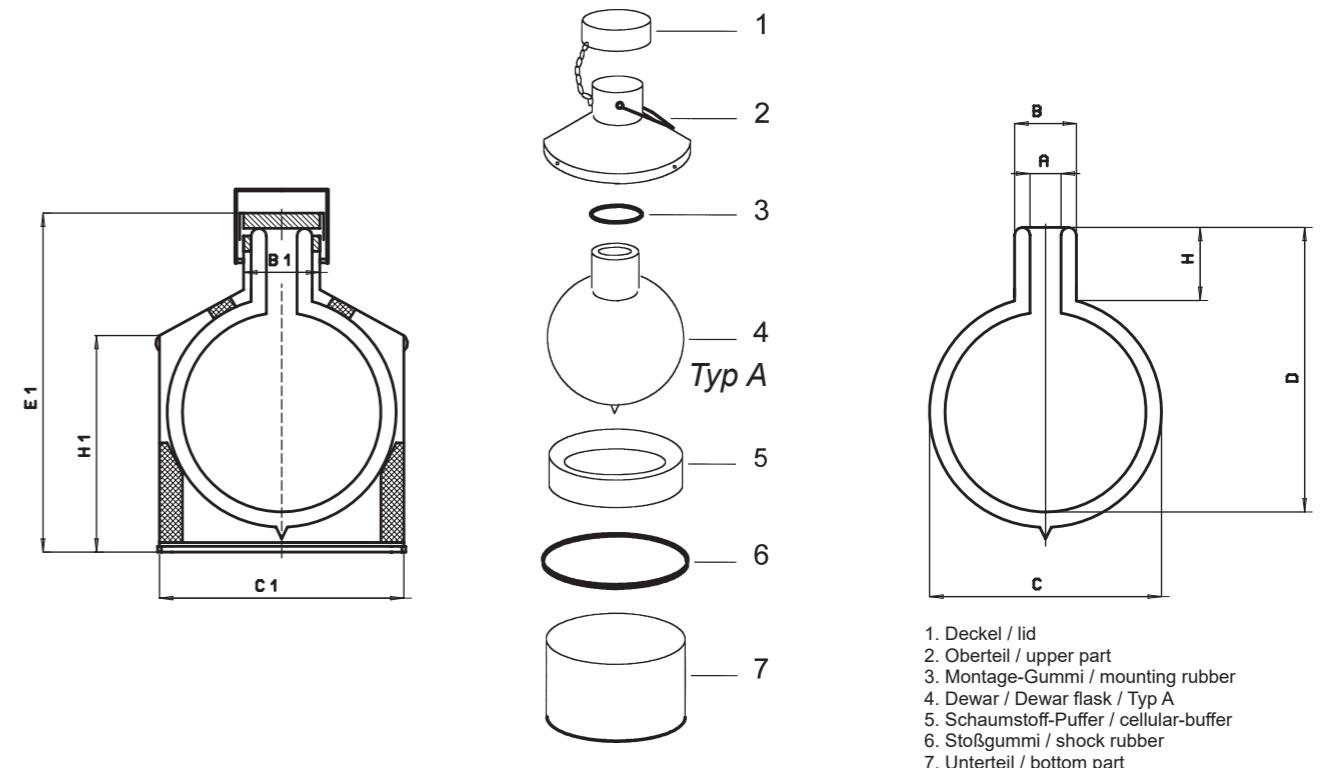
Typ T / Type T

Application

- Maintaining temperature of liquids
- The cooling agent heated by a thermostat is cirulating in the chamber of this double-walled vessel.

Advantages

- made of Borosilicatglass 3.3
- operating temperature max. 300 °C / Typ T max. 120 °C / Typ T-GL
- coolant ram pressure max. 0,5 bar
- high thermal resistance
- very good chemical resistance
- transparency
- deliverable:
- various special jacket joints
- special dimensions

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
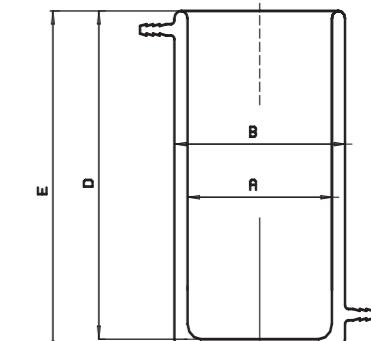
Different sizes and measurements on request.

Temperierbecher aus Edelstahl Tempering beakers made of stainless steel Bécher de temérisation de d'ancier fin

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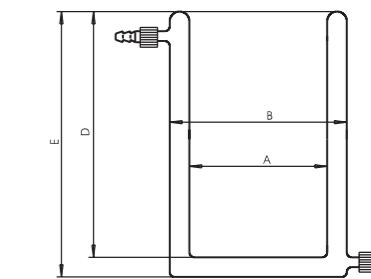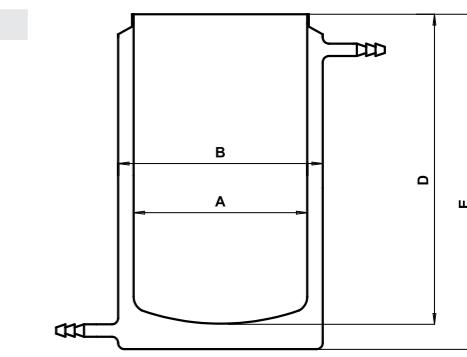
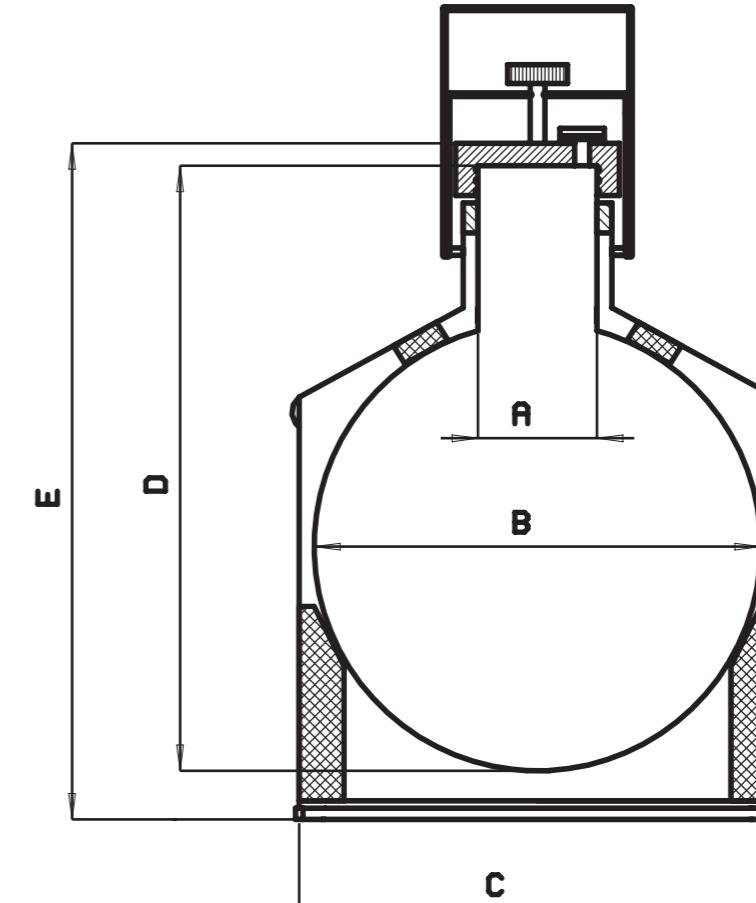
A user's manual and further designs can be found under
www.kgw-isotherm.com/produkte/kanne.html

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Applications

- storage and transport of aggressive waste in laboratory

Advantages

- good tipping stability
- made of Borosilicate glass 3.3 (DURAN)
- The glass is resistant against nearly all kinds of chemicals
- According DIN12111, 12116, 52322 the glass fulfils the conditions of:
 - hydrologic class 1
 - acid class 1
 - alkaline class 2
- The glass surface will only be attacked by:
 - hydrofluoric acid
 - heated phosphoric acid
 - alkaline solutions (in high concentration)

Standard equipment

- screw lock
- overpressure safety device
- handle
- glass insert
- red-coated cover
- warning stickers

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Application

- flexible applications in laboratory
- cooling of small samples
- storage and transport of :
- liquid nitrogen (LN2)
- dry ice (CO_2)
- other coolants
- suitable as insulating flask

Advantages

- high mechanical breaking strength
- easy handling
- made of stainless steel
- with carrying handle (except type DSS/GSS 500)

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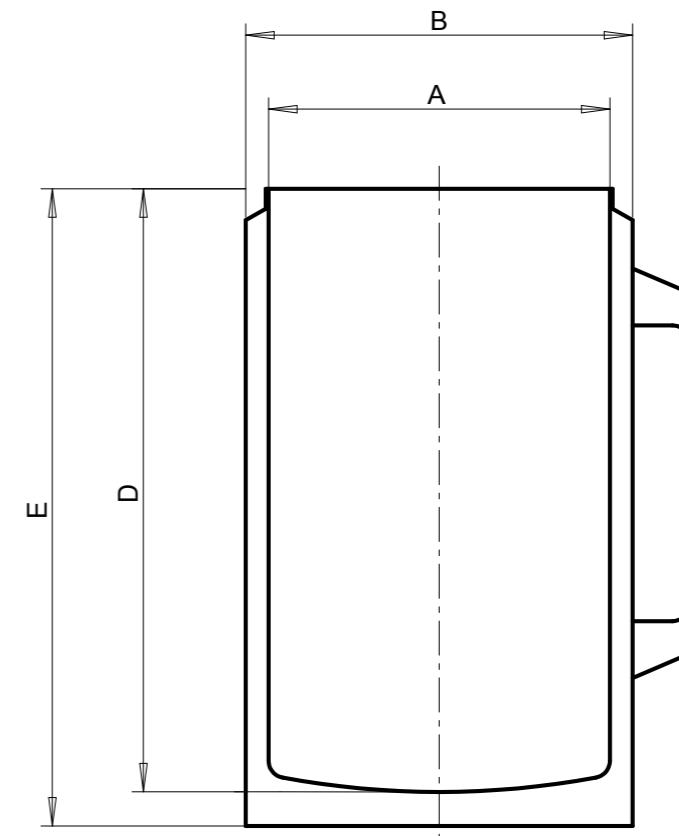
Thermal Insulation of Data Measuring Logs

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In refrigeration plants, conveyor, paint or temper ovens, or anywhere where extreme temperatures are encountered, the exact control of temperature is critical in the quality of the product. Warm handling processes are frequently assisted by temperature data devices which control the product through the high and low temperature boundaries. Usually the operating temperature of data logs is in the band of 0° to 60° Celsius. So a thermal insulation is necessary.

ITET is a series of powerful performance hot and cold protection devices developed to provide extreme control in critical situations. It is possible to control the parameters at the critical spots in the product while the process is running. The sensor is fixed at or in the product and measures the temperature. A remote cable is no longer necessary.

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* The operating time is the maximal time while the inner temperature is in the range of 0°C to 60°C ($T_{Start}=20^\circ\text{C}$). The operating time depends on the surrounding temperature

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Special Request

If you have any insulation problem, please send us your requests. We will try to offer an insulation device, that meets your requests.
For further information contact us for assistance.
(Germany) Fax. 0049 721/95897-77

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Applications

- gas cleaning by condensation e.g. in vacuum systems

Advantages

- made of Borosilicate glass 3.3
- good chemical resistance
- good thermal resistance
- deliverable:
 - Type GKF: aluminium foot (Standard), all around protection
 - with special joints: e.g. KF flange 16-25

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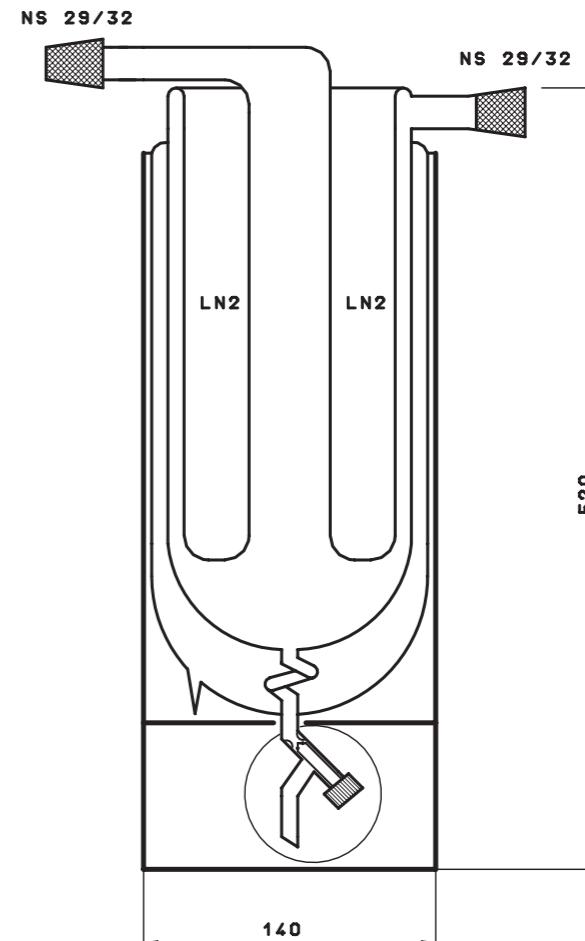
Construction Details

- two condensation walls, cooled by liquid nitrogen (LN2)
- integrated Dewar flask
- viewing stripes for observing the condensation process
- condensate outlet with O - ring sealing
- protective cover made of aluminium
- coolant capacity : ca. 1100 ml
- condensate capacity: ca. 250 ml
- joints:
 - to pump: ground cone NS 29
 - to plant: socket NS 29
 - insulating lid

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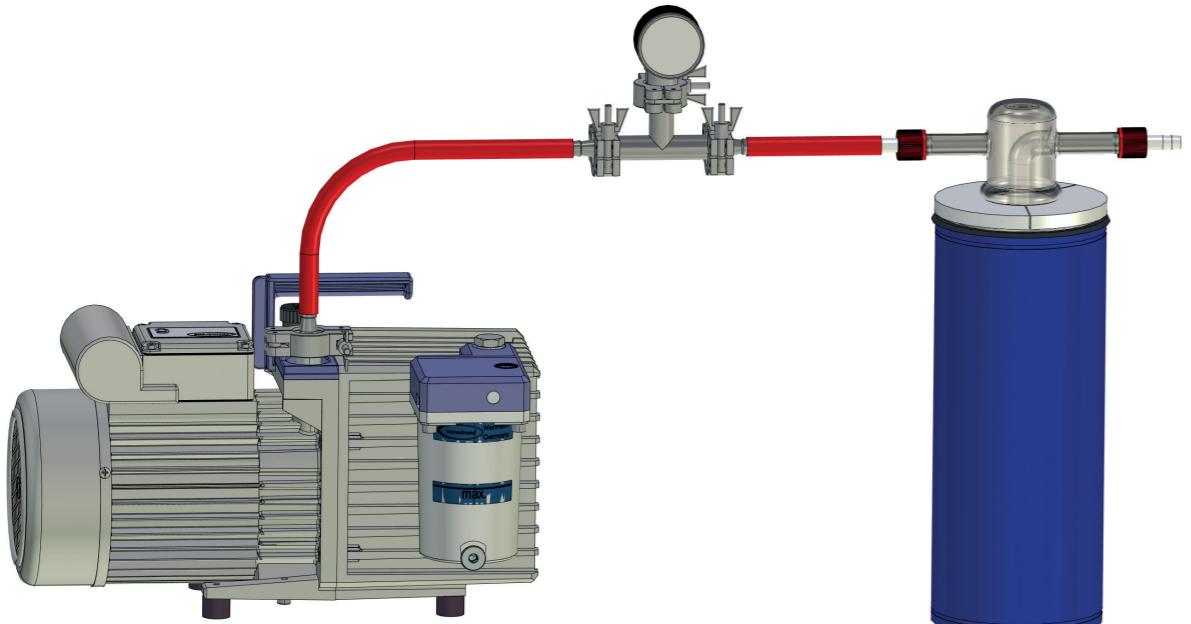
Applications

- gas cleaning by condensation
e.g. in vacuum systems

Advantages

- made of Borosilicatglass 3.3
- good chemical resistance
- good thermal resistance
- simple low cost design
- quick exchange is possible

deliverable:


OK = spherical joint S29 with O-ring sealing
K = spherical joint S29
GL = Glass screwthread with olive GL 18
- special dimensions
- special joints
e.g. small flange KF NW
DN flange

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CO₂ wire basket for cold traps


If dry ice (CO₂) is used as coolant, it might get difficult to exchange the glass cold trap. The CO₂ can fill up the existing space inside the Dewar flask. It is nearly impossible to replace the trap into the Dewar flask, then. Therefore we designed a wire basket for easily placing the cold trap into the already filled Dewar flask.

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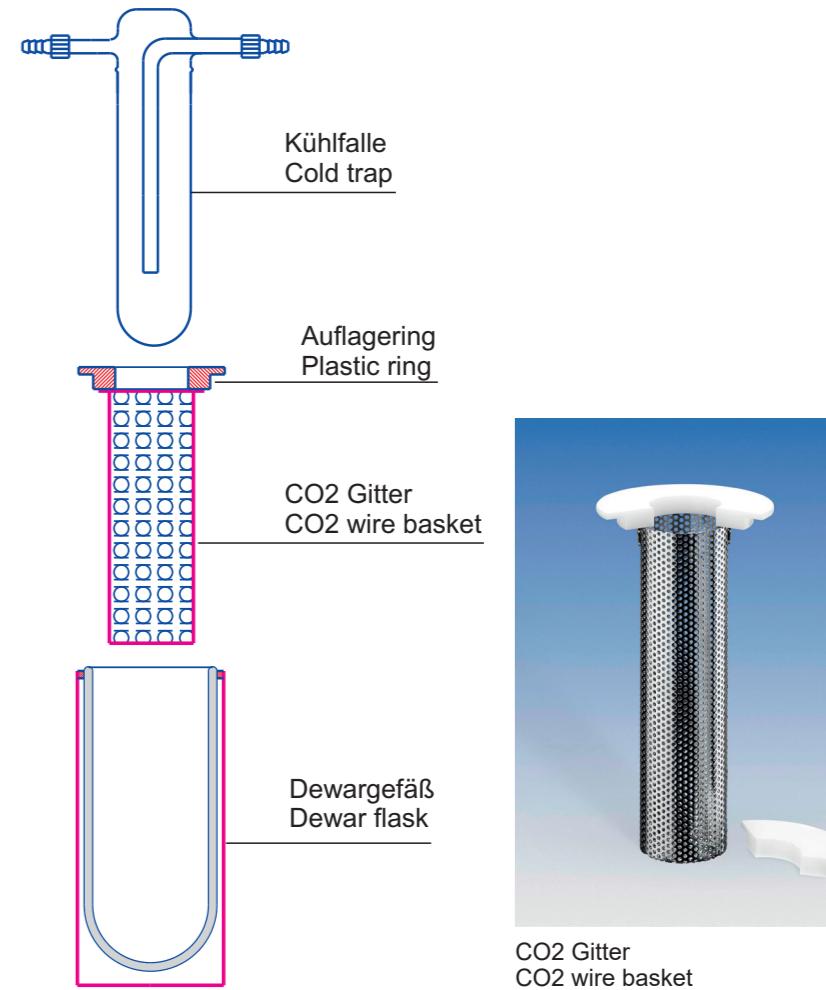
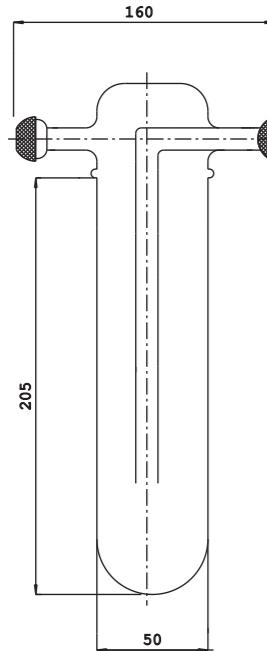
Kühlfalle zu KF 29-K
Cold trap for KF 29-K

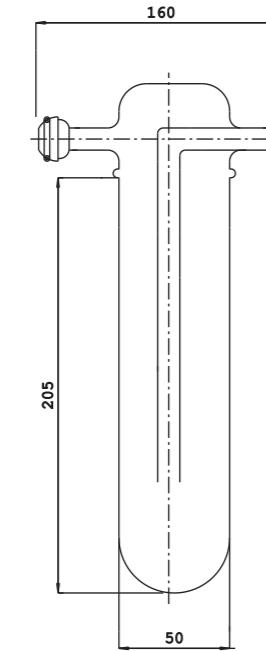Kühlfalle zu KF 29-OK
Cold trap for KF 29-OK

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
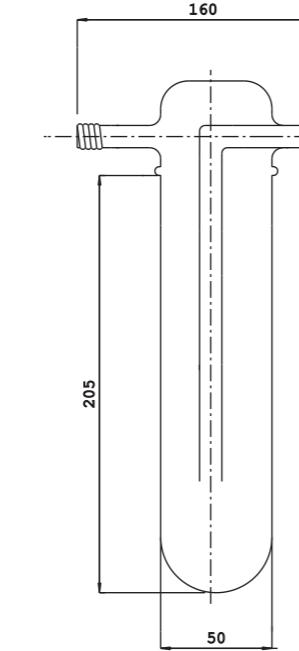
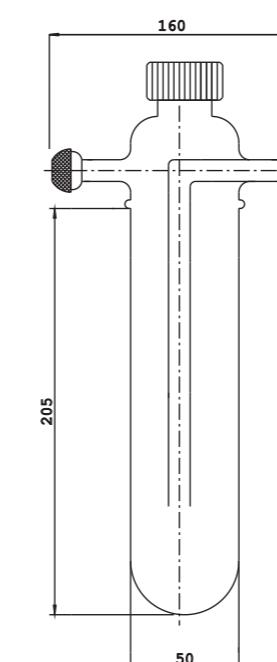
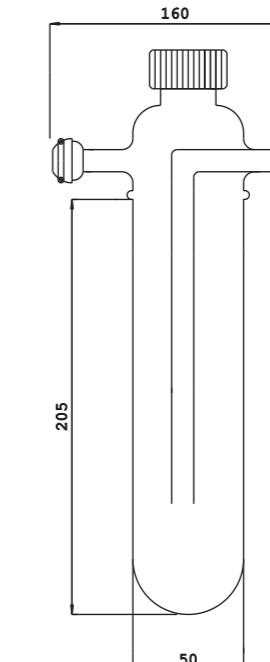
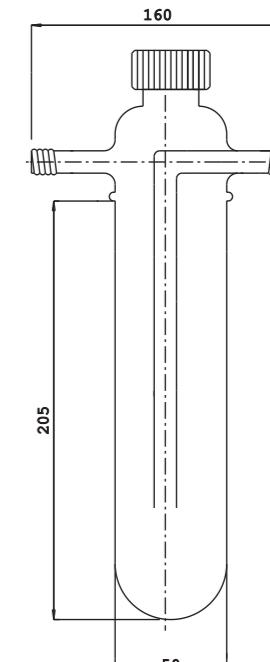
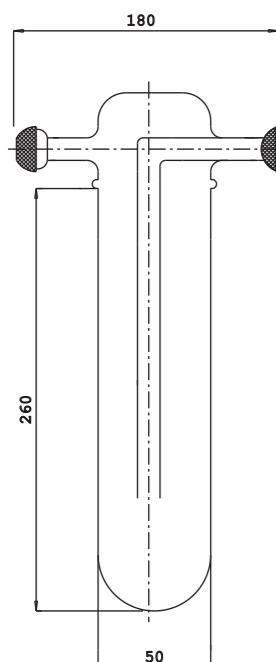
Kühlfalle zu KFL 29-K
Cold trap for KFL 29-K

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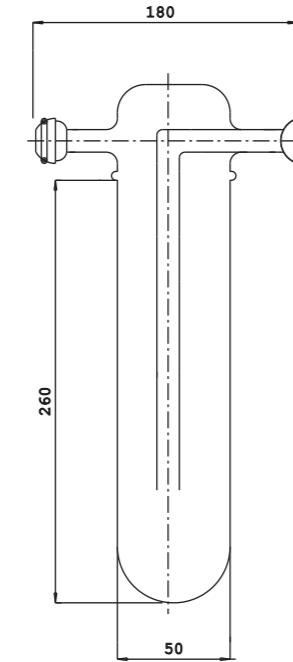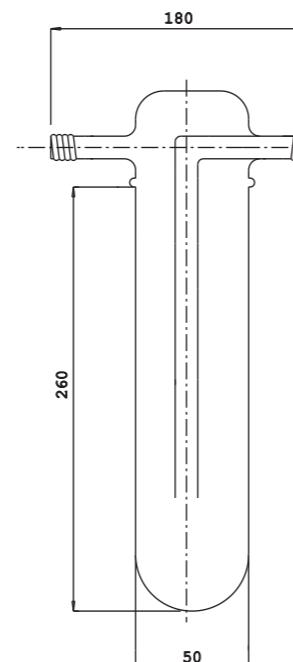
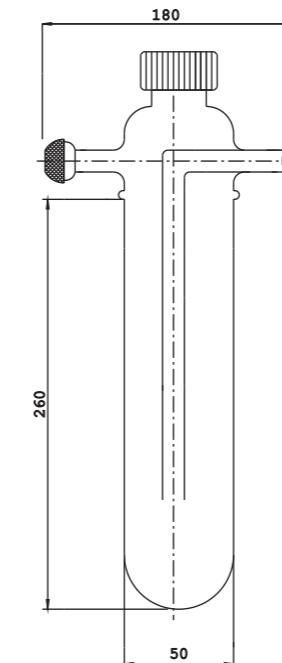
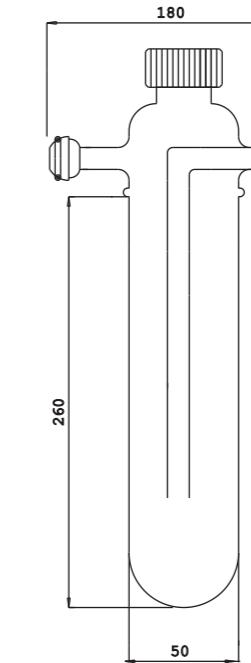
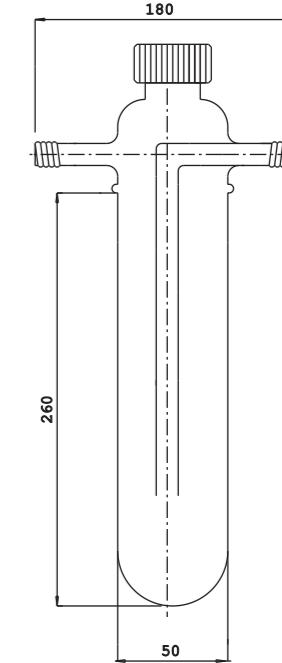
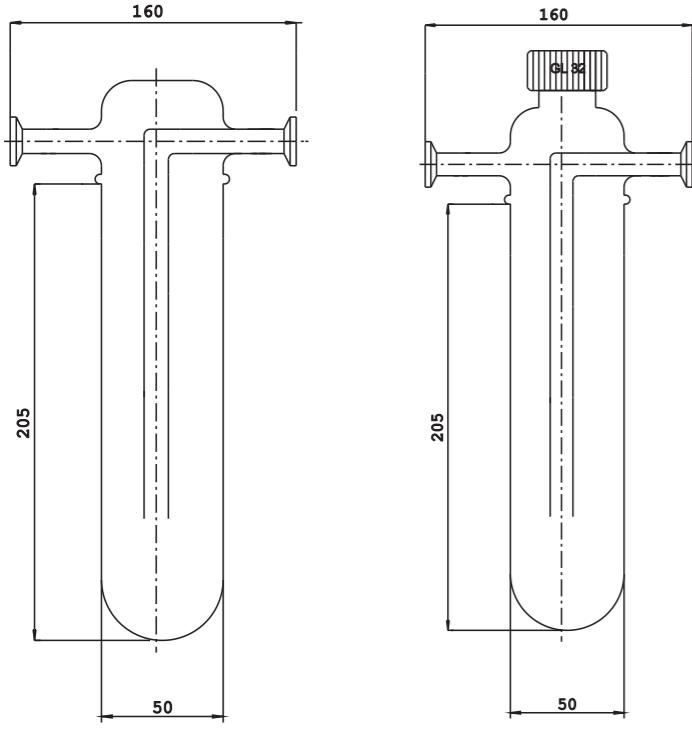
theoret.Vol. = Theoretisches Volumen

theoret. cap. = theoretical capacity

Cold Traps (cold finger condenser with KF NW flange)

Simple cold traps, which are also known as cold finger condensers, are used in standard applications in conjunction with vacuum pumps and the condensation is forced out of the humidity or solvents used in the chemical applications. The cold trap is used to protect the vacuum pump in this case and the speciality of this type of construction is that the user does not need to use a stand to hold the cold trap. The cold trap has a rim, which is hung inside the Dewar flask's support ring. The cold trap can be changed very quickly without any problems arising. With glass vacuum flanges you are able to connect this cold trap directly to the metal flanges of a high vacuum pump.

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
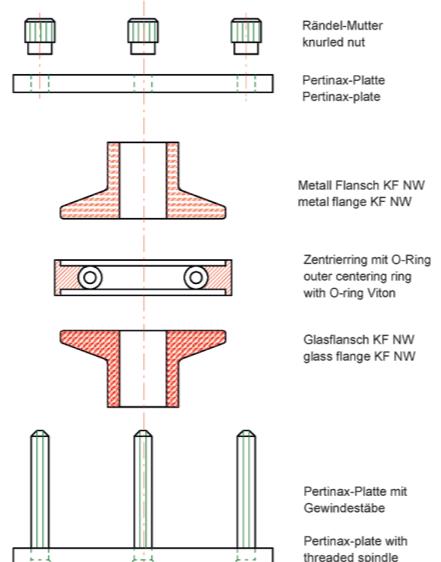
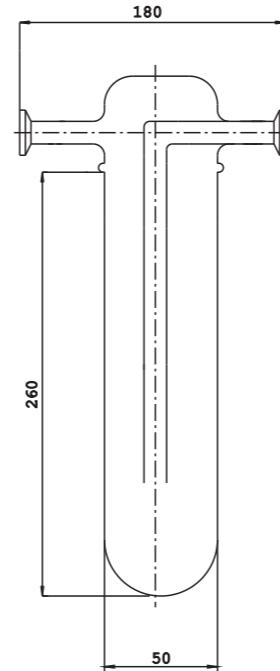
Kühlalfe für KFL 29-NW 16
Cold trap for KFL 29-NW 16

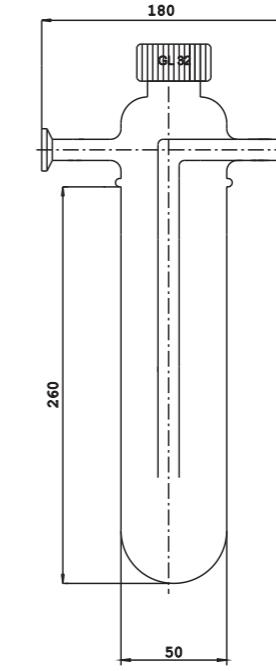Kühlalfe für KFL 29-NW 16-A
Cold trap for KFL 29-NW 16-A

Cold Traps long version (cold finger condenser with KF NW flange)

Simple cold traps, which are also known as cold finger condensers, are used in standard applications in conjunction with vacuum pumps and the condensation is forced out of the humidity or solvents used in the chemical applications. The cold trap is used to protect the vacuum pump in this case and the speciality of this type of construction is that the user does not need to use a stand to hold the cold trap. The cold trap has a rim, which is hung inside the Dewar flask's support ring. The cold trap can be changed very quickly without any problems arising. With glass vacuum flanges you are able to connect this cold trap directly to the metal flanges of a high vacuum pump.

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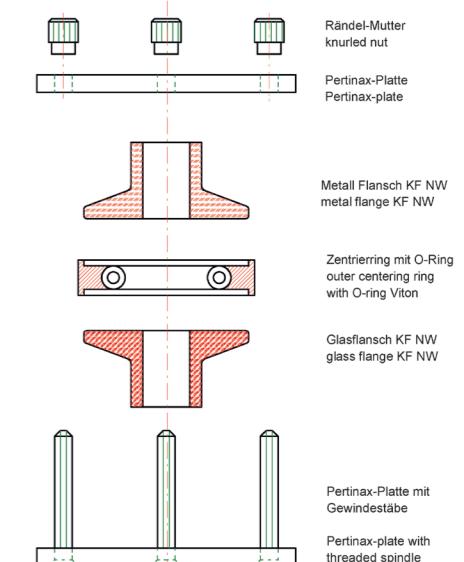Typ ISO-K100-Z **type ISO-K100-Z**

Characteristics

- Cold trap out of V2A stainless steel with a condensate capacity of approximately 1,5L
- baffles for optimized condensation
- reliable and easy handling
- no stand material for holding the cold trap necessary
- Dewar flasks made of glass according to DIN EN ISO16496 and a capacity of 16L, the stainless steel Dewar flask has a capacity of 18L
- protective casing of glass Dewar flask made of blue coated metal, aluminum stucco or stainless steel
- for liquid cooling agents, e.g. LN₂ (approx. -196°C)
- for solid cooling agents CO₂ with solvent (CO₂-wire basket necessary)

Article description

Dewar flask type 33C with suitable supporting ring and high performance cold trap S-K100-Z

Optional

Also available with CO₂ wire basket and Dewar flask 33CAL (with Stucco-aluminium) casing


Dewar flask out of stainless steel type DSS-D250/450 with suitable supporting ring and high performance cold trap S-K100-Z

Also available with CO₂ wire basket for dry ice cooling

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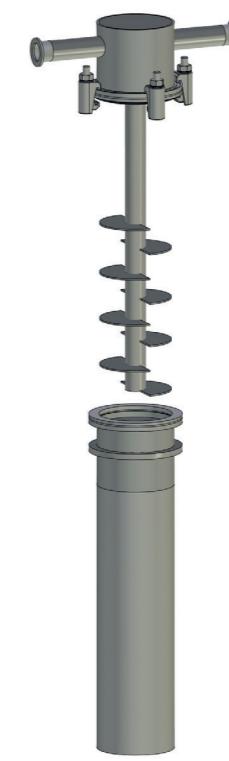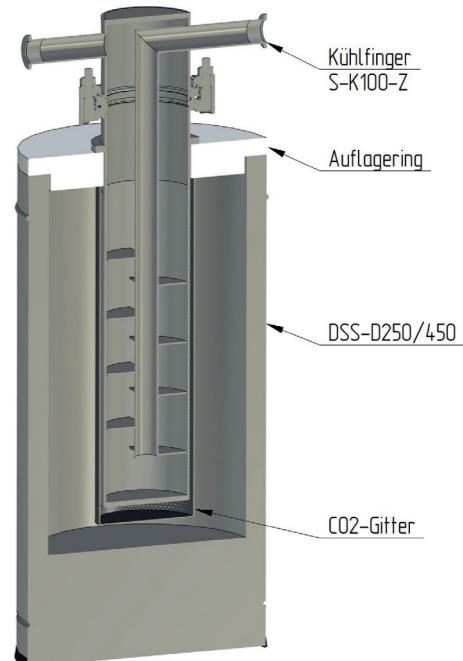
Typ ohne Manometer: CP 2
 Artikel Nr.: 1708

Applications

- flexible applications in laboratory e.g. for:
 - drying processes
 - condensation processes
 - vacuum technology
 - gas cleaning

Advantages

- compact design
- easy handling
- suitable for aggressive mediums
- glass components in Borosilicate 3.3
- O-ring sealed spherical joints incl. fork clamps
- mobile rack
- PE table board
- parallel operation with 2 cold traps. On failure of one cold trap the other cold trap keeps on working
- alternate operation with 1 cold trap. The cold trap in use is chosen by the valve adjustment

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
Typ GP3 mit Manometer
 Best. Nr. 17530

Chemistry pump device GP3 with cold trap and pump fork with aeration


Type GP3 is a mobile chemistry pump stand with a cold trap and a pump fork with a separate aeration cell. The pump fork has four independently switched vacuum valves and additionally four independently switched aeration valves, so that every sample holder can get individually evacuated or aerated. Through this type of assembling, every sample holder can get disconnected or aerated separately without affecting the vacuum of the other sample holders. Type GP3 has a small flange KF-NW 16 with vacuum manometer. Instead of the vacuum manometer you can also connect an electronic vacuum gauge. The pump stand will be delivered as standard without suitable vacuum pump.

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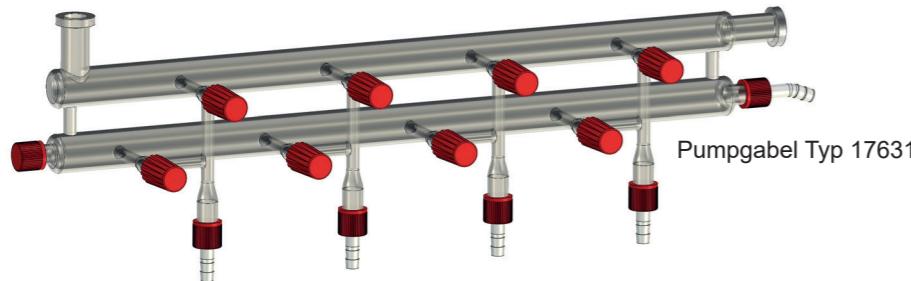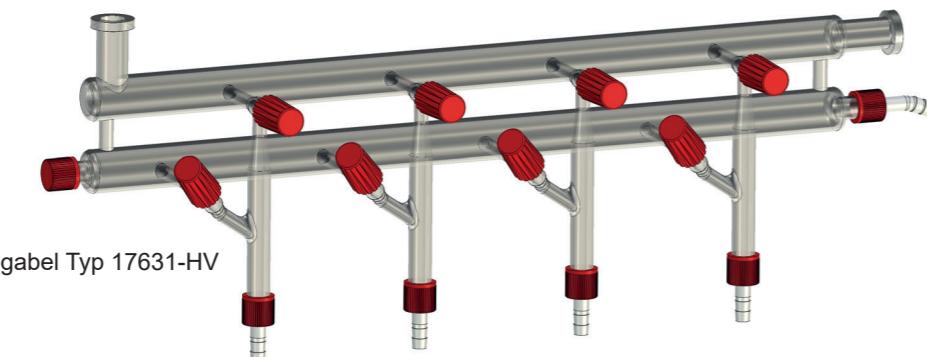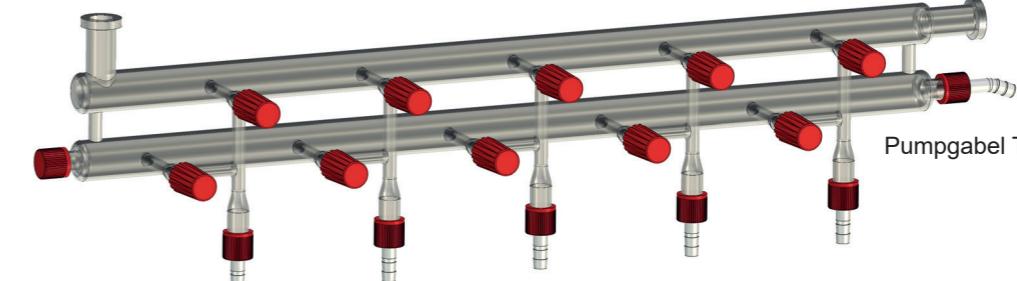
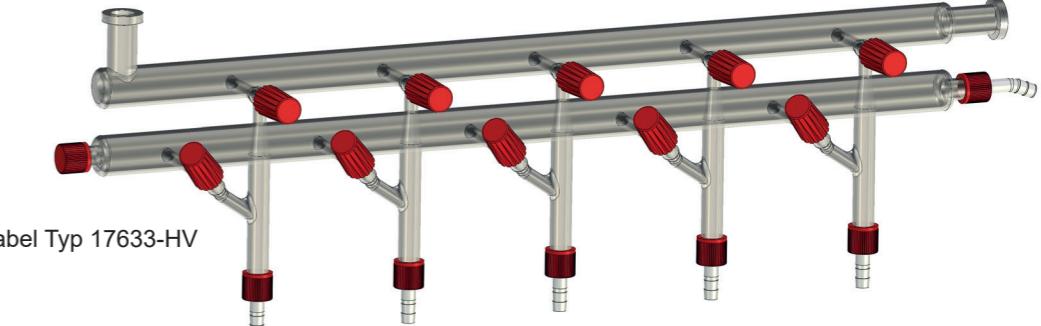
Area of application

For condensating water, solvents or gases in connection with a vacuum pump

Characteristics


- Reliable and easy handling
- Borosilicate glass 3.3 according to ISO 3585 (DURAN)
- Produram valve with o-ring valve
- Easy disassembly and cleaning of the valves
- Interchangeable valve stamp
- Individually configurable

Vacuum connections for the pump fork / schlenkline

- glass threads with screw-on cap and PTFE olive, e.g. GL18
- vacuum glass flanges, e.g. KF NW 16
- plane flanges, e.g. DN 25
- conical ground joints NS 29/32
- spherical ground joints S 29
- spherical ground joints Rotulex S 29 with O-ring seal

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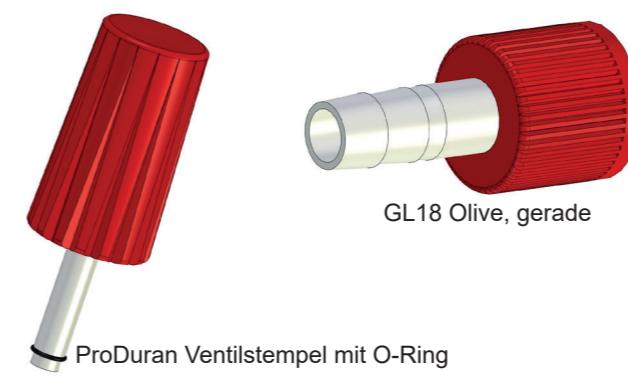Automatic level control of liquid nitrogen in a Dewar flask

In many cryogenic applications, such as shrinking of shafts and bushings in LN2, freezing of components or biological samples, using of LN2 as a coolant in cold traps, calibrating of sensors in LN2, etc., it is important to maintain a constant LN2 liquid level. This is achieved with the LN2 level control of KGW-ISOTHERM. With help of our LN2 level control, it is possible to keep a constant level of LN2 between the minimum and maximum sensor inside the flask. The KGW LN2 level control can be connected to almost every LN2 container.

The LN2 level control operates as follows:


If the LN2 level is sinking below the minimum sensor, a signal will be sent to the level control system. The LN2 level control opens the 24V LN2 magnetic valve, which is located at the transfer siphon which is connected at the LN2 storage container. Because of the over pressure inside of the LN2 storage container the LN2 will be pushed out the container and passes through a transfer line to the working Dewar Flask. The LN2 will float the working Dewar Flask until the maximum sensor will be dipped into the LN2 level. The maximum sensor is cooled as a result and gives a signal to the level controller. This will switch off the magnetic valve power supply so that the valve closes until the LN2 level will sink below the minimum sensor. Then the process described above will be repeated again.

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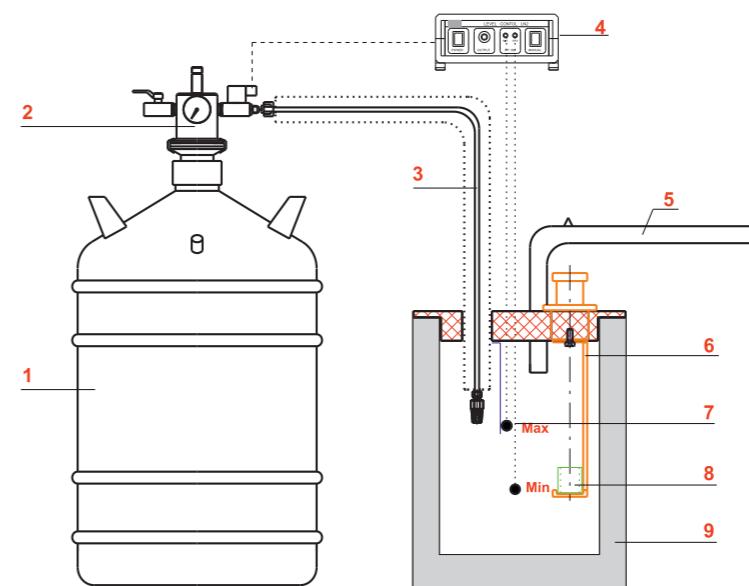Automatic LN2 - level control devices for manual shrinking

To make manual or automatic connections safely through shrinking, it is important that a sufficient cooling is ensured at the shrinking components. LN2 is usually used as coolant for shrinking bushes or shafts. With help of an automatic LN2 level control system the necessary cooling can be ensured. The cooling via LN2 can be secured for the shrinking process.

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
Automatic LN2 - level control devices for manual shrinking with two holder for the shrinking bushes, Dewar flask, lid, gas outlet tube and LN2 store vessel with 35 litre. Art.No. 2760-35

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
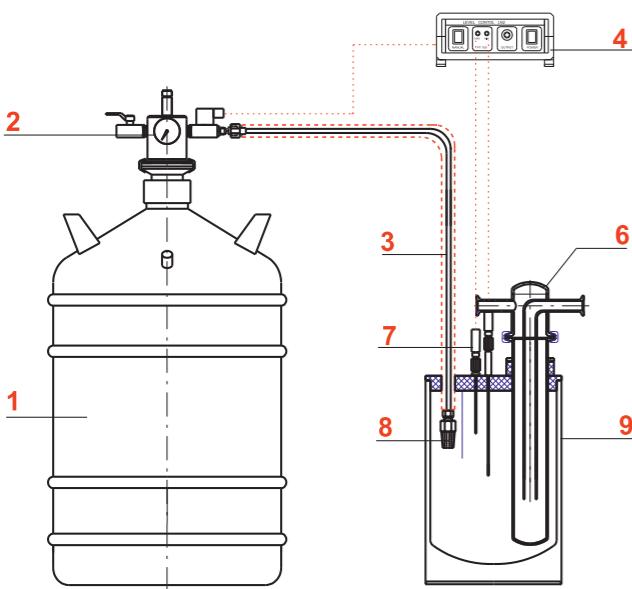
Automatic LN2 - level control with a stainless steel cold trap S 54V-K16-Z, Dewar flask, lid, gas outlet tube and LN2 store vessel with 25 litre. Art.No. 2750-25

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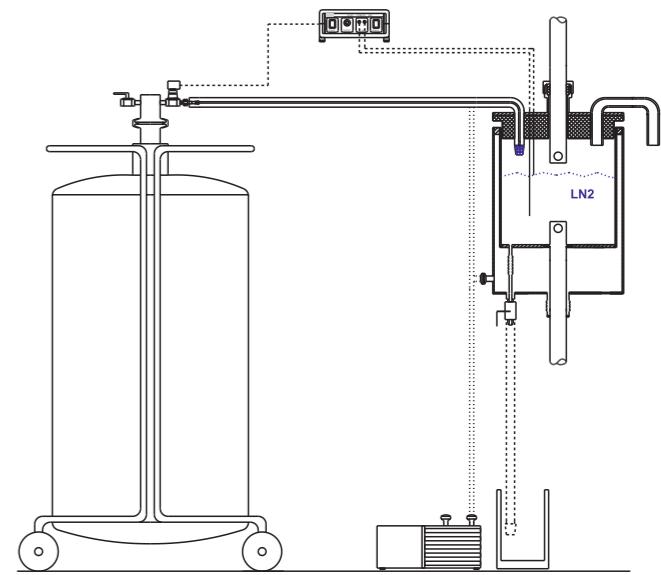Automatic LN2 - level control devices on a Dewar flask for a tension test

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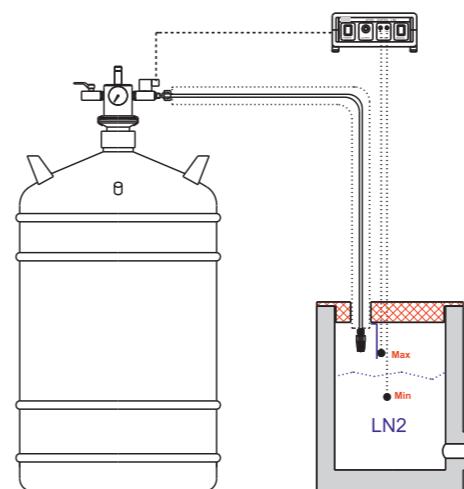
Automatic LN2 - level control devices with a corresponding Dewar flask

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Automatic LN2 - level control devices with two Dewar flasks

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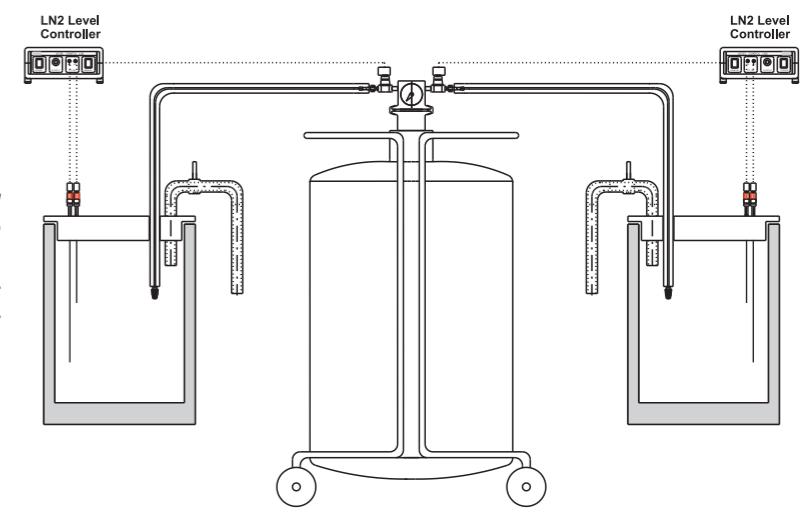Typ: LEVEL CONTROL LN2
Art.Nr.: 1301

Application

- Level control of liquid nitrogen (LN2) in unsealed Dewar flasks
- Automatic refill device for liquid nitrogen (LN2)
- For refilling, the controller opens the magnetic valve at the reservoir tank

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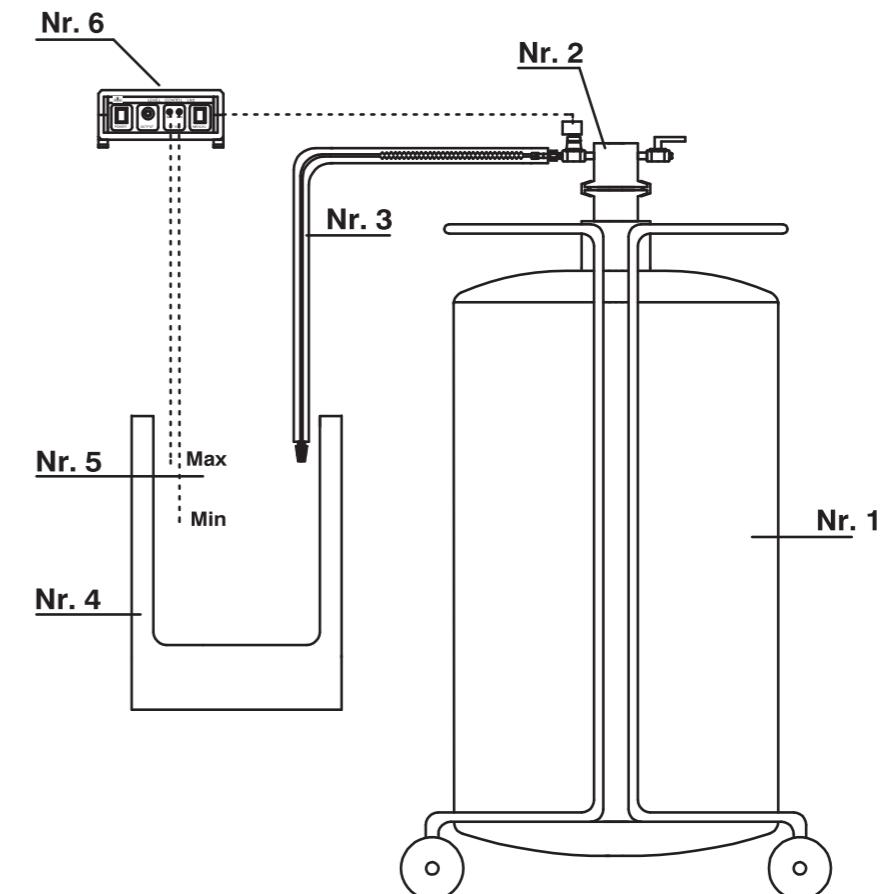
Accessories

- Temperature sensor PT 100
Art.Nr: 1302
- Magnetic valve for liquid nitrogen (LN2) (24 V)
Art.Nr: 1303
- LN2 transfer siphon with magnetic valve for reservoir tanks with safety device TYPE 1
Art.Nr: 1304
- LN2 transfer siphon with magnetic valve for reservoir tanks without safety device TYPE 2
Art.Nr: 1305

technische Daten siehe rechte Seite / technical data see right page / caractéristiques techniques voir page de droite

Automatisches LN2 Niveauregelsystem mit vier Regelführer Automatic LN2 level control system with four control sensors

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Application

- Level control of liquid nitrogen (LN2) in open or closed dewar vessels
- automatic refilling device for liquid nitrogen (LN2) with level monitoring
- For filling, the device opens and closes an LN2 supply solenoid valve on a storage and supply tank

Advantages

- manual or automatic filling via Min-Max-sensor
- with upper and lower alarm sensor
- cable breakage protection of the sensors
- optical indication of the switching status
- optical level indication by the control sensors
- acoustic and optical error indication
- one sensor lance with all 4 control sensors and one reference sensor

Scope of delivery CRYO LC




- Control unit with connection for LN2 solenoid valve (230V AC)
- option with connection for LN2 solenoid valve (24V DC) and sensor function monitoring

Article number:

Typ/e: CRYO LC for 230V
LN2-magnetic valve,
Art. Nr./ Art. No.:4711-230

Typ/e: CRYO LC for 24V
LN2-magnetic valve,
Art. Nr./ Art. No.:4711-24

Sensor lance, position of the control sensors and length of the lance according to customer specification
Art. No 4711-L

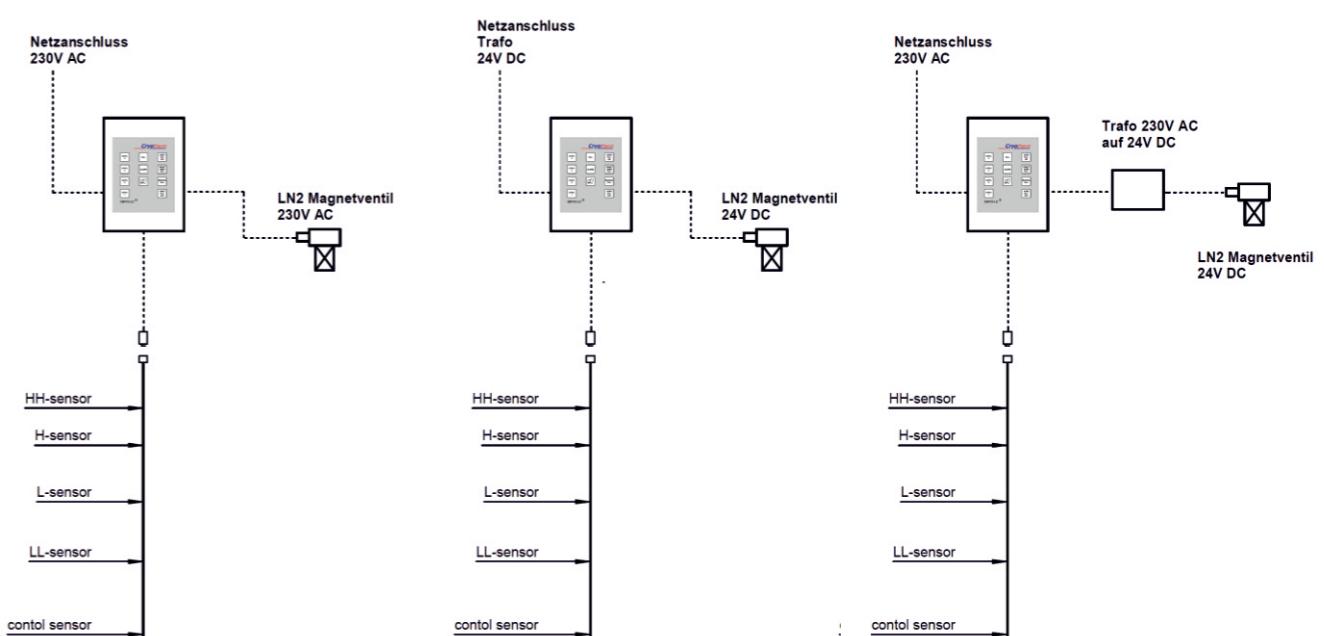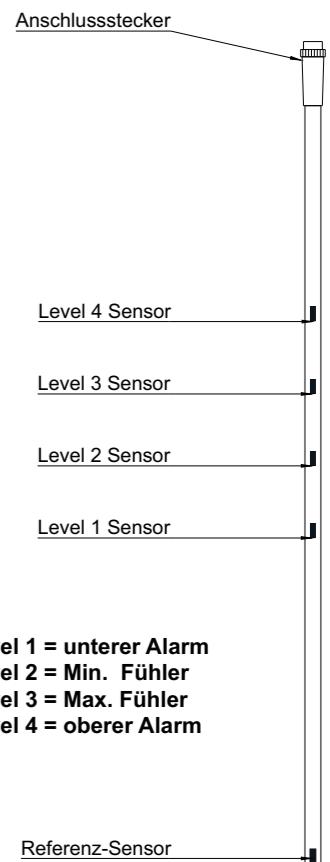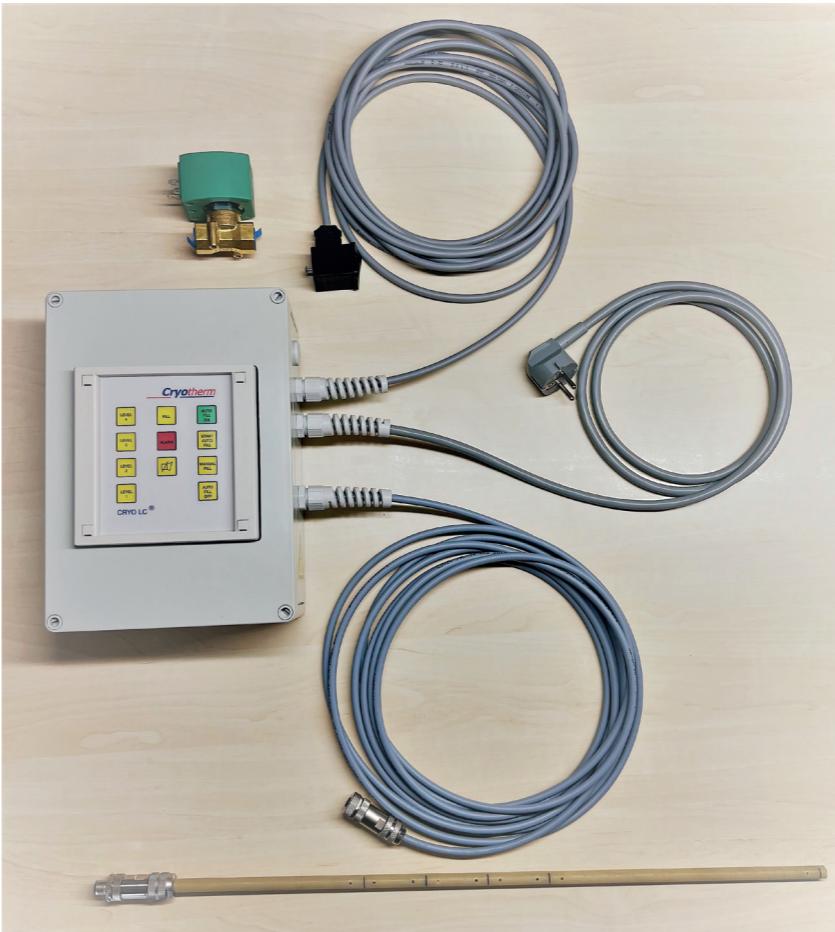technische Daten siehe rechte Seite / technical data see right page / caractéristiques techniques voir page de droite

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
Filling system

The filling system is an equipment to fill in LN2 into a vessel without removing the LN2 siphon.
The filling system is made of:

- 1) The filling flange with LN2 valve, waste gas valve, screw adapter, O-ring and clamp for NW 50.
Filling flange **No.: 2606**
- 2) The filling hose has got two screw couplings for the screw adapter.
Filling hose **No.: 2607**

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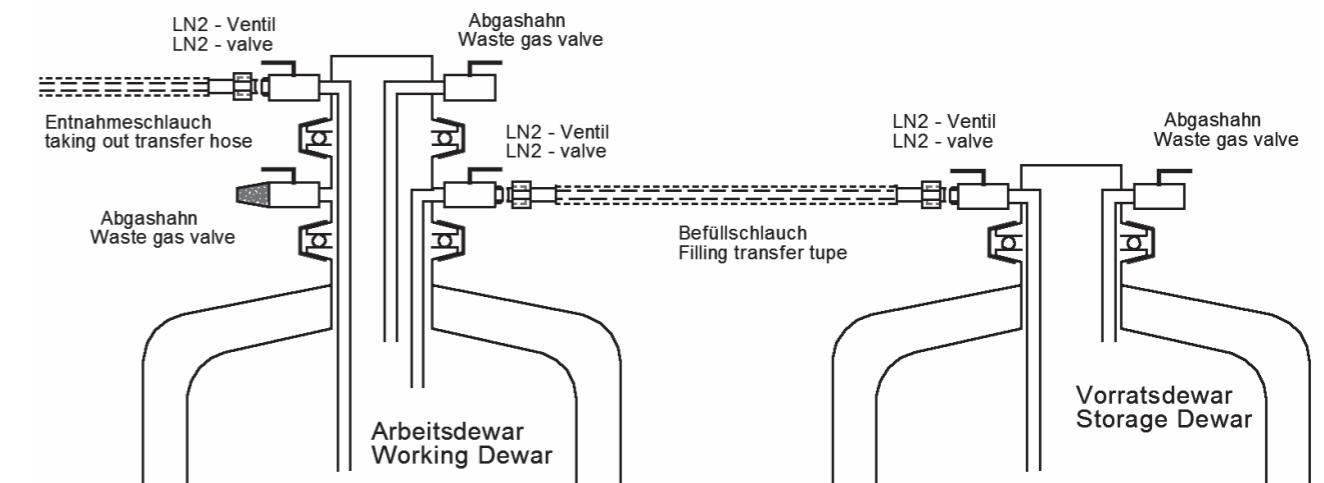Characteristics

High mechanical stability, the container is made of stainless steel. Low evaporation rate by multi-layer superinsulation in the vacuum space. Very durable by robust constructive design, production and choice of material. Long-lasting insulating vacuum by the use of getter materials.

Standard equipment

Integrated safety devices at the container neck
Vacuum lock with safety valve
Easy running castors
Small-flange connection KF NW 50
Transfer siphon with transfer hose covered with metallic tissue
Hand and protection ring
Contents gauge digital
Pressure attachment by fluid removal without external power supply

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easy
robust
efficient
low-priced

Juno 30
with KF NW 50 flange,
over pressure valve
and manometer

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Area of application

- for liquid nitrogen (LN2)
- pressure-less space for coolants
- laboratory technology
- medicinal technology
- biotechnology
- storing and transport

Features of performance

- reliably and economically
- low own weight
- compact construction
- nominal capacity 30 litres

Types and accessories

- double wall, vacuum insulated aluminium vessel
- multi-layer super insulation in the vacuum space
- carrier ring with side grips
- loose lid-on plug
- vessel with KF NW 50 for a LN2-siphon
- over pressure valve
- manometer

Temperature range
LN2 -196°C

Pressure range
Pressure less or up to max. 1,3 bar

Safety advises and regulations

- always wear protective glasses and protective gloves
- national regulations for laboratories
- company- guide lines
- safety regulations for handling with liquid gases
- ADR regulations

Evaporation rate
Daily evaporation and static holding time were determined at 20°C, 1013 mbar, container stationary, lid closed and without transfer siphon. These are nominal values that may change depending on the age and use of the container.

Flüssigstickstoff-Behälter aus Edelstahl Typ JUNO, (30 L) Liquid nitrogen container made of stainless steel Type JUNO, (30 L)

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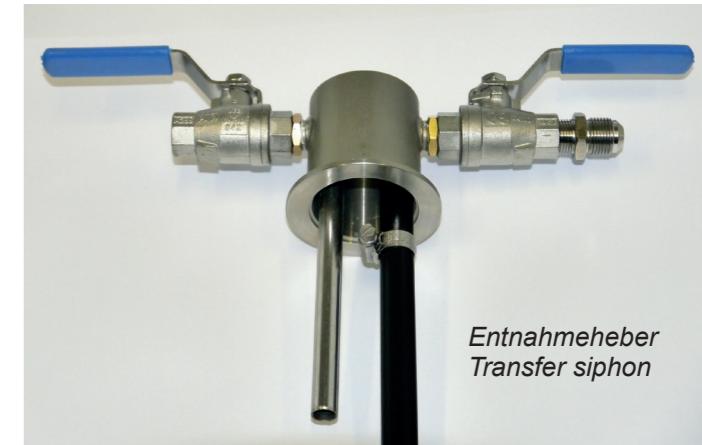
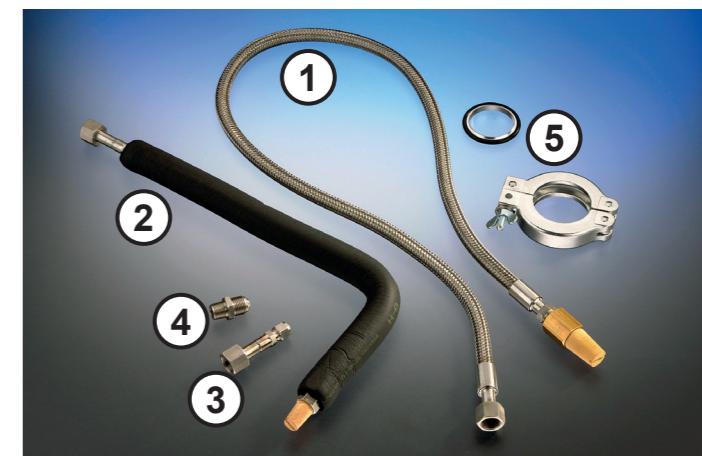
Cylindrical Dewar vessels made of stainless steel for liquid nitrogen for use in laboratories and industries

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- + Manufactured out of high quality stainless steel
- + Surface with a brushed finish
- + Use of super-insulating foil and insulating lid
- + Easy carrying due to side grips
- + High stability
- + Flattened internal bottom


Dewar vessels Type DSS - D with insulating lid and side grips
Dewar vessels Type DSS - DF with flange on request
Intermediate sizes and special designs can be supplied on request

Dewar vessels Type DSS - D with insulating lid and side grips
Dewar vessels Type DSS - DF with flange on request
Intermediate sizes and special designs can be supplied on request


*Zylindrische Dewargefäße aus Edelstahl
Cylindrical Dewar vessels made of stainless steel*

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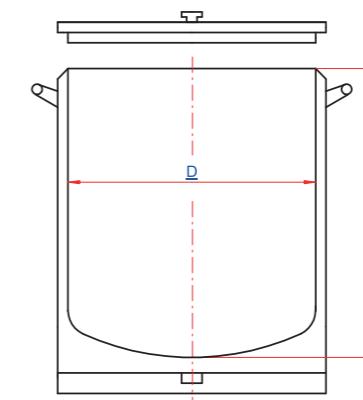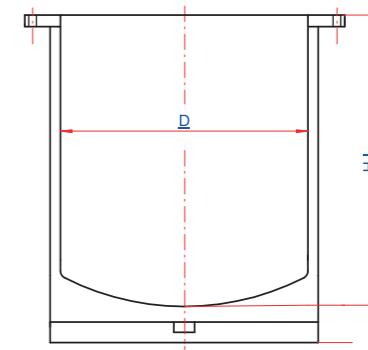Characteristics

Light aluminium container with compound materials in the neck range
Low own weight
Super-insulation foil in the vacuum space
Low evaporation rate

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Standard equipment

Handle or side grips
Vacuum valve
Small-flange KF NW 50
Loosely resting plug

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
Transfer siphon incl. safety valve, manometer, pressure reduction valve made of stainless steel, fluid valve made of stainless steel and screw adapter
Best.-Nr. 2611

① Entnahmeschlauch aus Edelstahl / Transfer hose made of stainless steel
Best.-Nr. 2612

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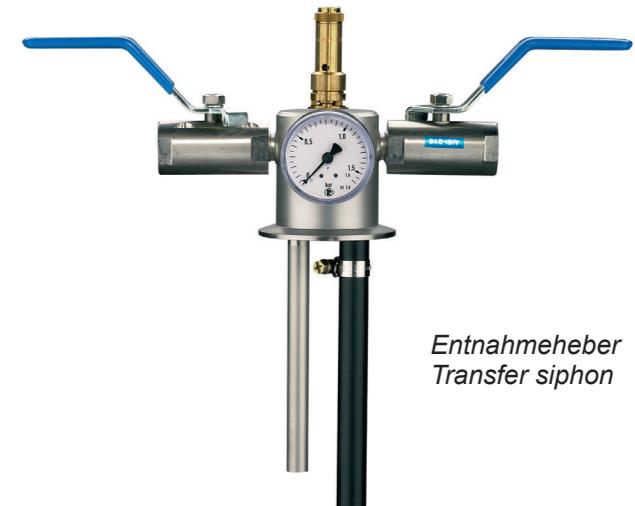
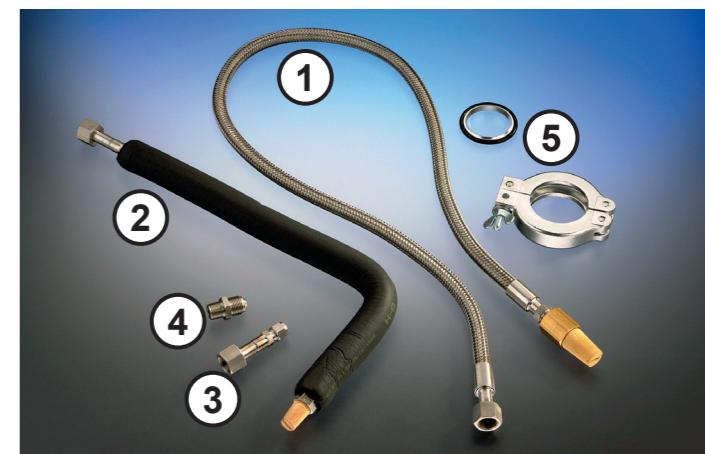
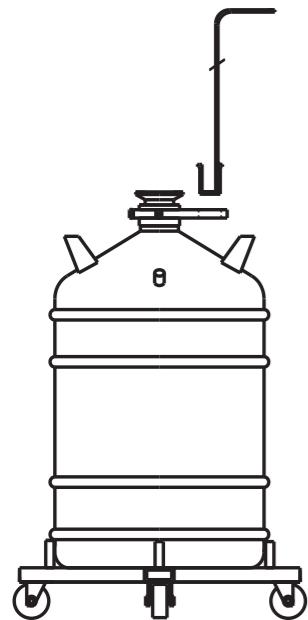
Safety advises and regulations

- always wear protective glasses and protective gloves
- national regulations for laboratories
- company- internal regulations
- safety regulations for handling with liquid gases

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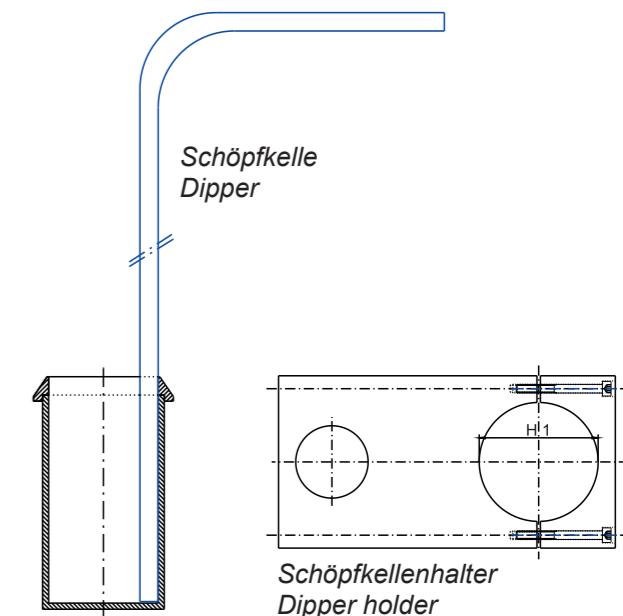Area of application

- pressure-less space for coolants
- Laboratory technique
- medicinal technique
- Biotechnology
- storing and transport

Features of performance

- reliably and economically
- low own weight
- simple LN2 removal
- easy cool down a sample in LN2
- nominal capacity from 12 to 35 litres

Description of Aluminium Dewar vessel Type DMT

Types and accessories

- double wall, vacuum insulated aluminium vessel
- multi-layer superinsulation in the vacuum space
- side grips or handle
- loosely resting plug
- vessel with KF NW 50 and dipper holder
- dipper to take out LN2

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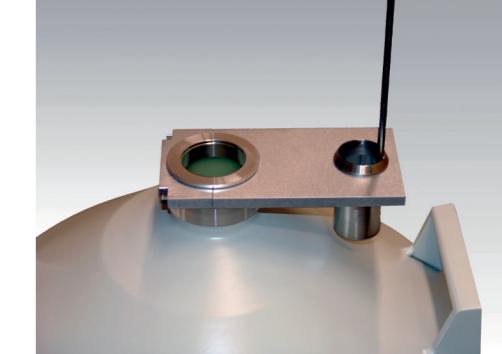
Characteristics

Light aluminium container with compound materials in the neck range
Low own weight
Super-insulation foil in the vacuum space
Low evaporation rate
All the vessels comply with international regulations applicable to the transport of hazardous goods by road (ADR) or rail (RID).

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Standard equipment

Transport handle
Vacuum valve
Loose lied-on plug
Inclusive canister
Canister made of stainless steel, standard
Canister made of plastic, on request

technische Daten siehe rechte Seite / technical data see right page / caractéristiques techniques voir page de droite

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Daily evaporation and static holding time are measured at 20°C, 1013 mb, vessel stationary, lid closed and without inventory system. These values are nominal and can vary according to age and application of the container.

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Characteristics

Light aluminium container with compound materials in the neck range

Low evaporation rate

Super-insulation foil in the vacuum space

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Characteristics
Light aluminium container with compound materials in the neck range
Low evaporation rate
Super-insulation foil in the vacuum space
Low own weight
The liquid nitrogen is absorbed in a porous substance
No risk of spilling liquid nitrogen
The samples are preserved in gaseous phase
Complying with international regulations applicable to the transport of hazardous goods by road (ADR), rail (RID) and by air (IATA-OACI)

Standard equipment
Safety vacuum valve
Inclusive canisters
Handle or carrying strap

Accessories
Protective casing for transport

technische Daten siehe rechte Seite / technical data see right page / caractéristiques techniques voir page de droite

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Daily evaporation and static holding time are measured at 20°C, 1013 mb, vessel stationary, lid closed and without inventory system.
These values are nominal and can vary according to age and application of the container.

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What is KALTGAS ?

KALTGAS is a tempering system for a temperature range in between -180° C and $+100^{\circ}\text{ C}$ with a separation of the primary tempering devices and the specimen that is to be cooled. Low temperatures are available by the use of cold nitrogen gas. The gas is generated by evaporation of liquid nitrogen, then tempered by a tempering module and guided through an insulated pipe to the specimen.

Due to the separation of the primary tempering device and the specimen chamber, it is possible to operate various test chambers with only one KALTGAS system. Another feature of KALTGAS is the cooling of specimens in the free space. With this, the specimen is blown in the free space. It is therefore possible to test components or units while assembled or in operation.

Examples of applications

- Thermal testing of plastics, metals, compound materials and so on.
- Cooling of electronical components
- Tempering at tests of the quality control
- Shock-freezing of biological specimens, food and other materials
- Tempering of test specimens at
 - + tension- torsion, or notched bar tests
 - + chemical or physical tests
 - + process engeneering
- Specimen-cooling at diffractometers, spectrometers und detectors
- low-temperature techniques and cryogenics

Design

KALTGAS systems work with liquid nitrogen as cooling agent and have therefore the advantage of a very high cooling performance. A KALTGAS system is designed by the following components:

- Control unit
- Liquid nitrogen container
- Nitrogen evaporator
- Insulated pipe
- in the need a tempering chamber
- Tempering module

Due to the modular concept of KALTGAS, every system can be designed exactly according to the costumer-related requirements.

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KALTGAS system type TG-LKF for a temperature range from $+100^{\circ}\text{C}$ to -180°C , with a temperature stability of $\pm 0,2^{\circ}\text{C}$.

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KALTGAS system type TG-LKF-H for a temperature range from $+100^{\circ}\text{C}$ to -180°C , with a temperature stability of $\pm 0,2^{\circ}\text{C}$.

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KALTGAS system type TG-RID 40/50 for a temperature range from $+100^{\circ}\text{C}$ to -160°C , with a temperature stability of $\pm 0,2^{\circ}\text{C}$.

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KALTGAS system type TGK-40/50 for a temperature range from $+80^{\circ}\text{C}$ to -160°C , with a temperature stability of $\pm 0,2^{\circ}\text{C}$. The chamber has got a big vision panel.

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Typ TG-LKF-H

**Kaltgas cooling system
for cooling applications from +100°C (+248°F) to -180°C (-292°F)**

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
KALTGAS is a tempering system that utilizes the very low temperature of liquid nitrogen as a cooling agent. The Jet vaporizes the liquid nitrogen in a cryogenic container, producing a steady cryogenic gas flow. This LN2 vaporizer (Jet) can be adjusted to minimize liquid nitrogen consumption and to change both the cooling rate and the volume of the gas flow. The cryogenic gas flow is then piped through a vacuum insulated flexible metal line (N2 gas line) to a heat exchanger (heater). The heat exchanger is designed to heat the cold gas flow to the desired temperature. The final product, a clearly quantified gas flow exiting the heat exchanger at a clearly defined temperature, is ready for use as a cooling agent. A KALTGAS system needs only a few minutes to produce a cold gas flow with a temperature of -170°C (-274°F). Configuring the KALTGAS system to any customer's individual environment is possible thanks to various temperature controllers that can be embedded inside the safety control unit. Those temperature controllers are available with or without interfaces. With the standard safety control unit (SC5), a temperature stability of better than $\pm 0.2^\circ\text{C}$ can be achieved. The steady gas flow can be used to temper a chamber (TG-K KALTGAS).

In addition to their high cooling speed and good control stability, another advantage of KALTGAS systems is their modular design. By swapping out individual modules such as the N2 gas line, the LN2 vaporizer (Jet) or the heater, it is possible to change the cooling speed, LN2 consumption as well as the application. The basic modules, including the LN2 container, the vacuum pump and the safety control unit, remain unchanged.

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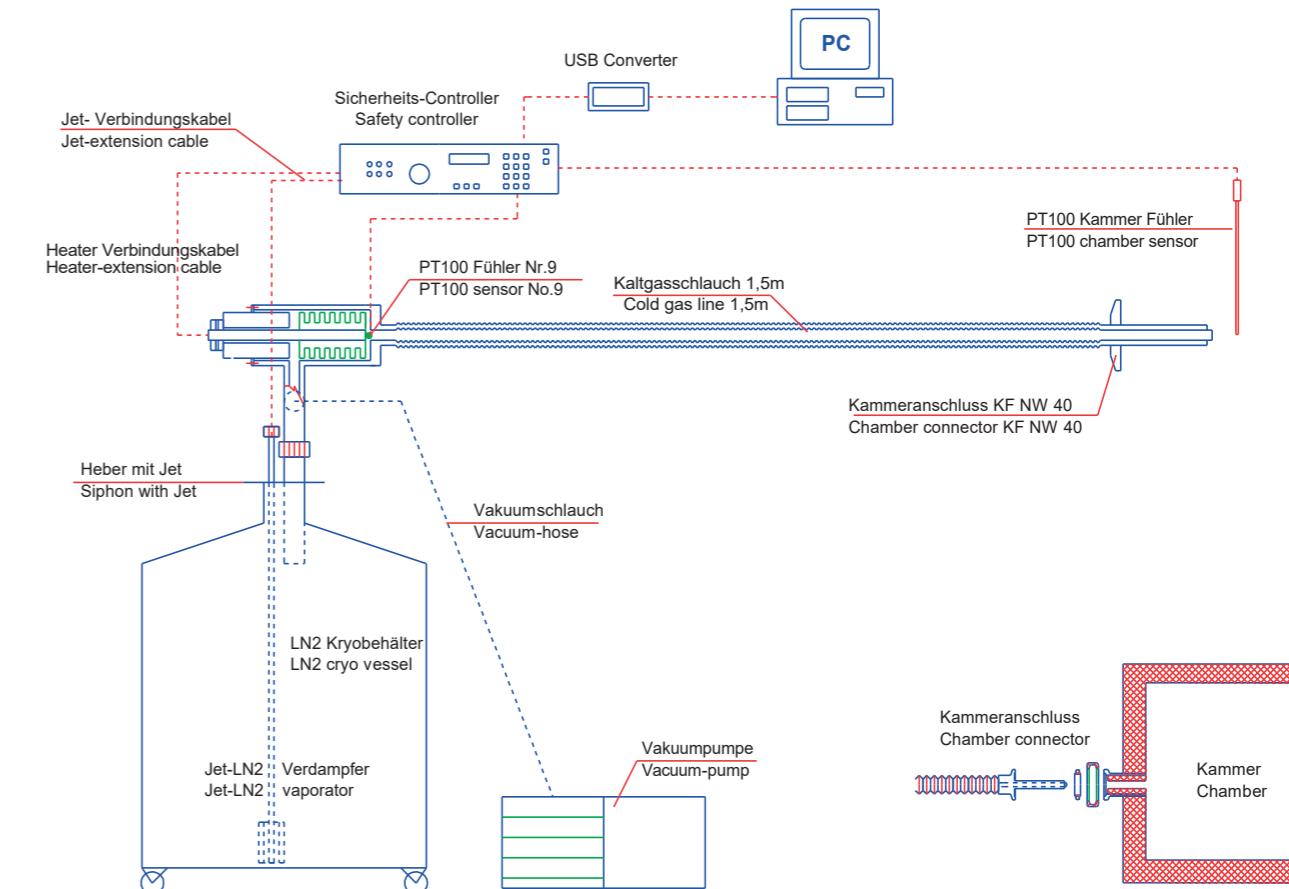Technical data:


Model TG-LKF-H 63/50

LN2 vaporizer (Jet)= 500 watts
Heater = 630 watts
LN2 consumption = 1.1l/h to 11l/h (liters per hour)
N2 gas line = V2A, length 1.8 meters, flexible,
with vacuum pump and gas outlet nozzle = KF NW 40,
siphon for LN2 container with KF NW 50

Order No.: TG-LKF-H 63/50

Accessories: LN2 container with 20 to
350 liter capacity

**Further information on KALTGAS systems available.
Technical data subject to change.**

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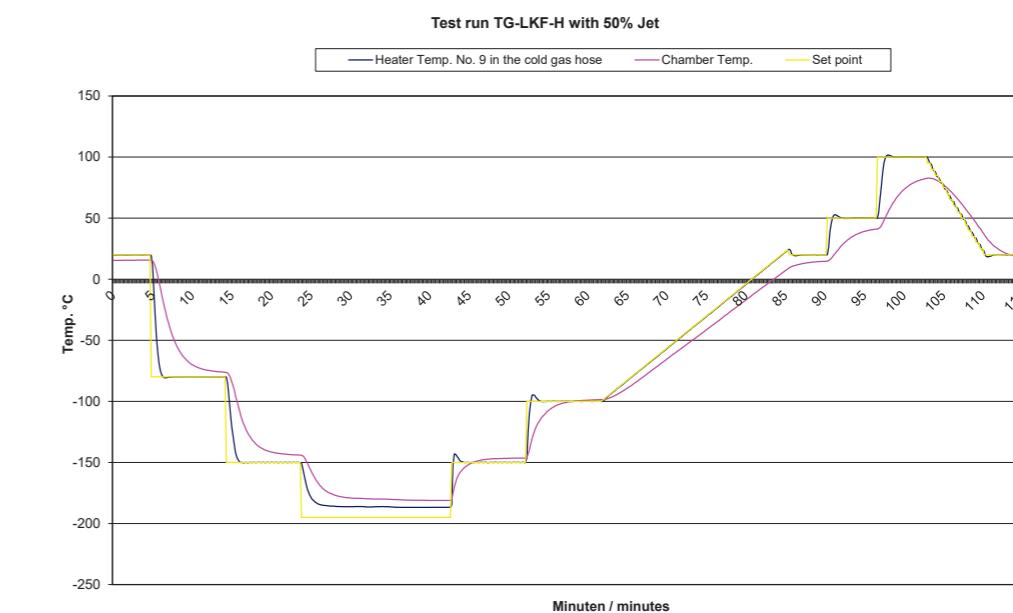All reaction vessels manufactured by KGW-ISOTHERM are based on plane flange system by Schott (NW 60 to NW 200). This standardised flange system provides customers with a comprehensive range of accessories and replacement parts. The O-ring groove in the vessel's flange facilitates simple and safe mounting of the vessel with the associated lid.

Our reaction vessels are offered according to standard in 3 versions: single-walled, double-walled and triple-walled. Vessels with, or without bottom outlet. Such bottom outlet can e.g. be a valve or a spherical joint.

There are different thermal-jacket connections possible for heatable reaction vessels. Beside the classic plane flanges DN15 or DN 25 also glass threads with plastic screw olives are offered.

All reaction vessels including the standard types are only made on customer's order. Therefore, modifications to joints or flanges are usually offered at no extra cost.

For reaction vessels with an indicated operating overpressure a temperature difference inside the glass wall of 50 K (for three-wall vessels 30 K) may not to be exceeded. For more information please have a look into our operating manuals under download on our website.

Each manufactured reaction vessel is compression tested after its production. Therefore the reactor must stand the one and a half time of nominal pressure at ambient temperature.

Other versions than those described here are always possible. On request we can supply you with a professional overview of our products per CAD drawings. A comprehensive range of accessories completes our program of reaction vessels. This is so large that it cannot be reproduced in the catalogue. If you need parts that are not shown, please do not hesitate to contact us. We will try to offer the desired part or a suitable alternative, then.

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
What do you need for a complete reactor with accessories

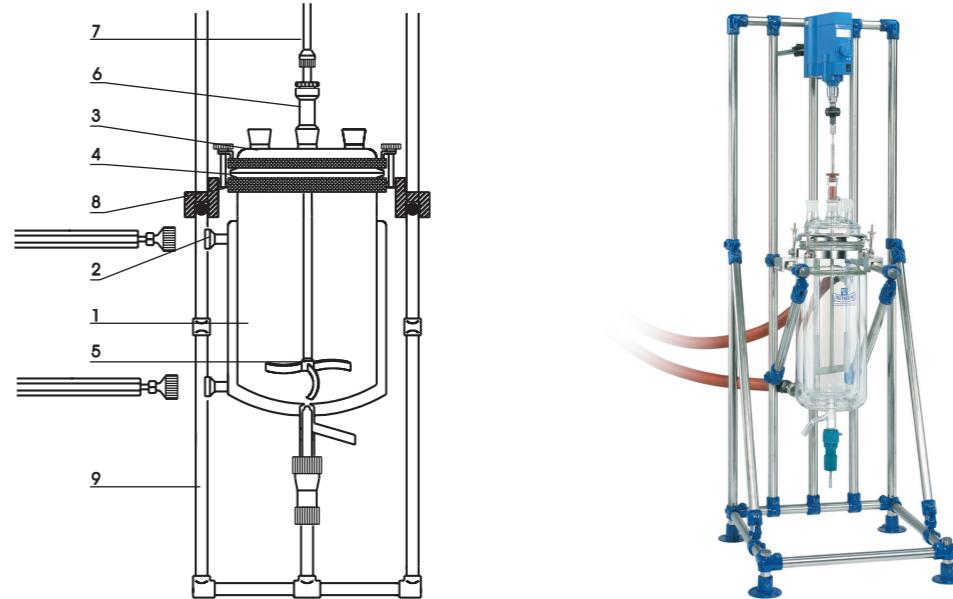- 1) **reactor**
(page 94 to 103)

capacity
size of flange NW 60 to 200
outlet valve (valve type L or type J)
thermal jacket connection (GL 18 or DN 15)
- 2) **reactor lid**
(page 104 and 105)

size of flange NW 60 to 200
centre socket (NS 29/32 or NS 45/40)
side socket (NS 14/23 to NS 45/40 / straight or angled)
- 3) **thermal jacket connection DN 15**
(page 107)

metal adapter DN 15
sealant DN 15 (silicone or PTFE)
thermal jacket hose (0,5 to 3 meters)
- 4) **O-ring**
(page 108)

size NW 60 to 200
silicone
FEP-coated
- 5) **seal and holdersize**
for reaction vessels
(page 108 and 109)

NW 60 to 200
type "Vario Grip"
type universal supporting device
type supporting device
- 6) **stirrer**
(page 110 and 111)

shaft of stirrer (8,10 or 16 mm)
anchor stirrer
propeller stirrer
- 7) **coupling for stirrer**
(page 112)

shaft of stirrer (8,10 or 16 mm)
type RKS
type RKG
- 8) **dynamic seals for stirrer**
(page 113)

joint (NS 29/32 or 45/40)
standard dynamic stirrer seal
PTFE-stirrer seal
magnetic stirrer coupling
type RKK
- 9) **frames**
(page 114 and 115)

table frame
stationary frame
movable frame cylindrically
Easy Frame

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
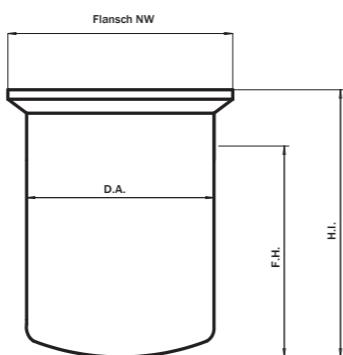
Reaction vessels cylindrically with "L" valve

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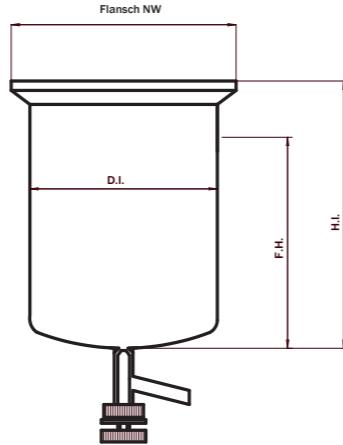
Reaction vessels cylindrically with 10 mm "J" valve

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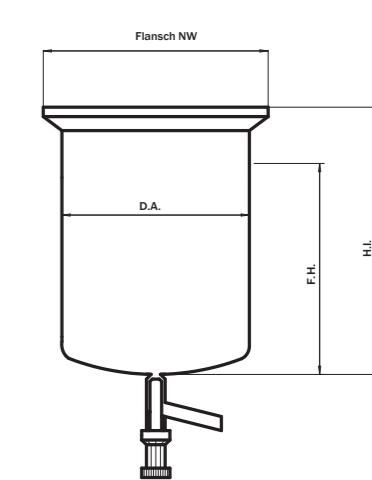
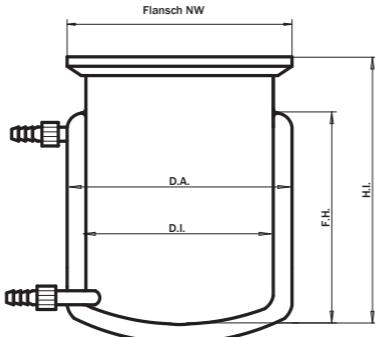
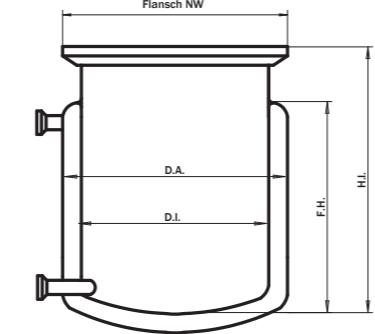
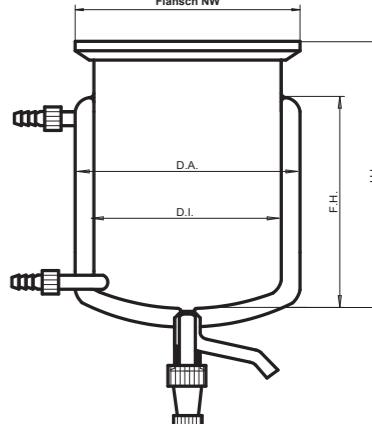
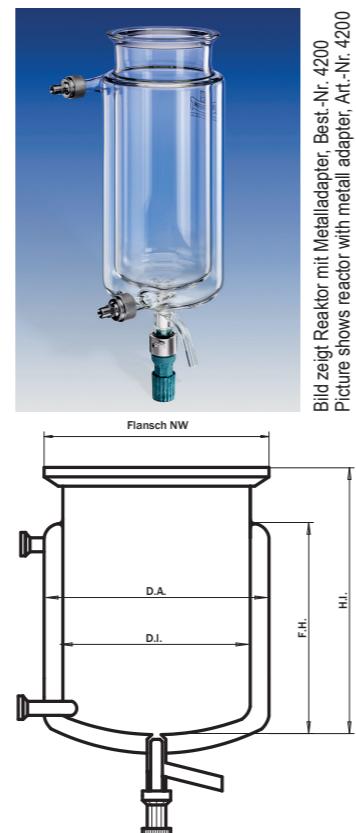
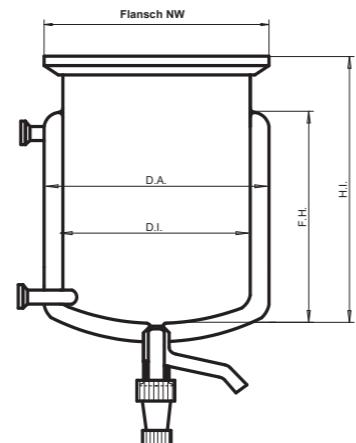
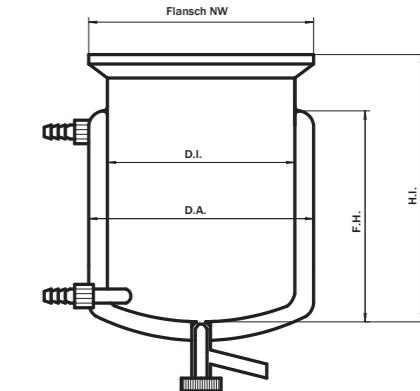
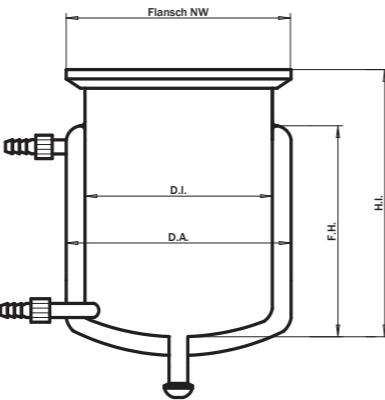
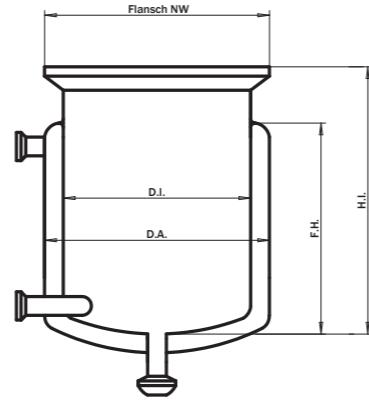
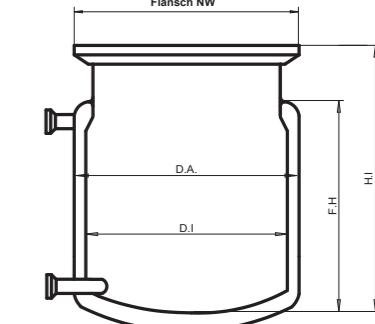
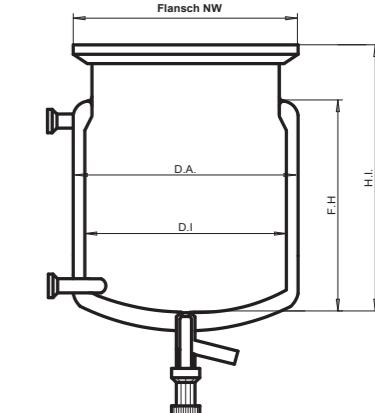
Heatable reaction vessels with flat bottom and thermal jacket joint GL18-olive

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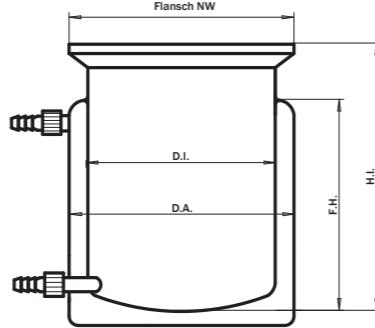
Heatable reaction vessels with flat bottom and thermal jacket joint DN15 flange

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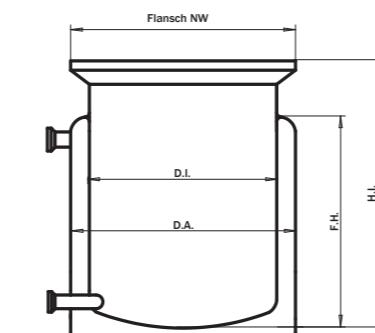Flansch und „J“ 10 mm Ventil


Heatable reaction vessels with vacuum jacket, with thermal jacket joint DN15 flange and "J" valve 10 mm

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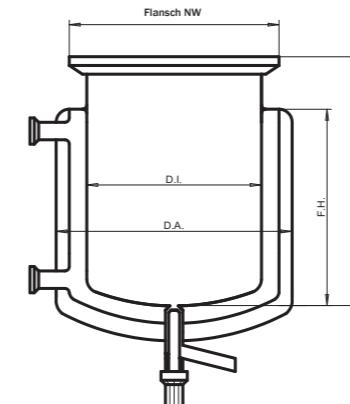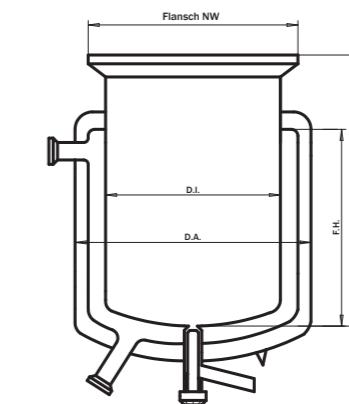
Heatable reaction vessels with vacuum jacket, with thermal jacket joint DN15 flange and 20mm „J“valve


Art.Nr	Vol.	D.A.	D.I.	F.H.	H.I.	Flansch	B.D.bar	Vers.

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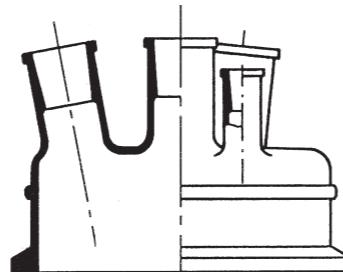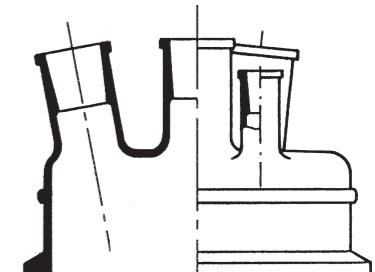
All lids are manufactured on customer's order, only. On request there can be alternatively set GL-glass threads, spherical joints or plane flanges instead of mentioned NS ports, partly without additional charges.

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Kugelschiff S 40

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Bottom outlet

The bottom outlet is available in three versions.

Industrial valve J (linear closing)

The industrial valve is a dead volume free PTFE spindle valve with sloped outlet. The spindle shall be pressed in without turning. In the valve housing is an additional spring mounted in. So the spring pressure of the spindle to the valve seat can be increased spring-loaded and the spindle regulates itself. This guarantees a good tightness at large variations of temperature. The industrial valve is available with a 10 mm bore hole (standard) and as desired with a 20 mm bore hole.

Laboratory valve L

The laboratory valve is a dead volume free PTFE spindle valve with the spindle being pressed in while turning. The valve has an 10 mm bore hole.

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Thermal jacket connection

The thermal jacket connection is available in three versions.

1 Glass-thread connection GL 18

- with plastic olive (PBT/PTFE) for hoses with an internal diameter of 10 mm.
- Internal diameter of olive is 7mm
- Accessories for thread connection GL18
- Plastic olive straight 4044
- Plastic olive angled 4044/A

2 DN15 glass flange

- Accessories for DN15 glass flange
- Metal adapter for thermostat hoses (DN15 to M16x1) 4200
 - Silicone sealant DN15 4201
 - PTFE-sealant DN15 4202

3 DN25 glass flange

- Accessories for DN25 glass flange
- Metal adapter for thermostat hoses (DN25 to M16x1) 4205
 - Silicone sealant DN25 4206
 - PTFE-sealant DN25 4207

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
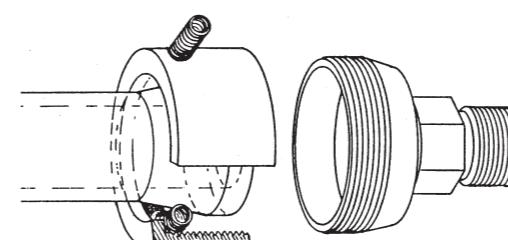
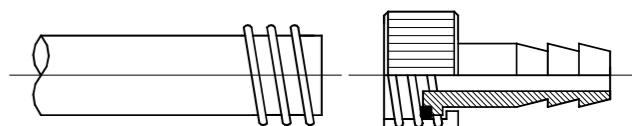
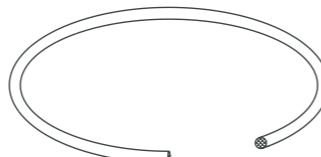
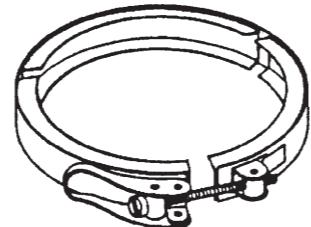
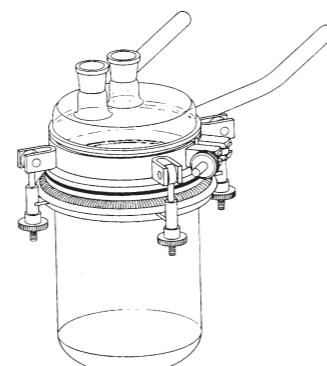
This holder allows to remove the lid or vessel from the holder or the set-up without dismantling the set-up. The supporting rods can be utilised vertically or horizontally.

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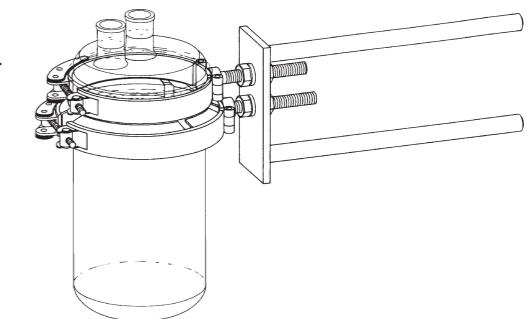
massive version for installation in a rod frame

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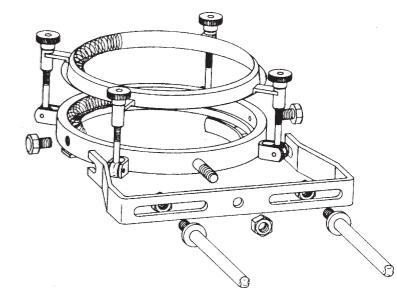
massive version for installation in supporting frame, thickness supporting frame : standard 27mm, on request 13 or 30mm

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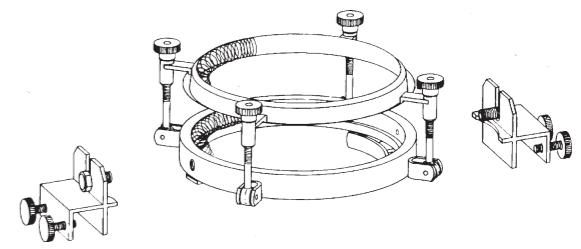
Following specifications are needed (please have a look at general drawing and example at bottom):

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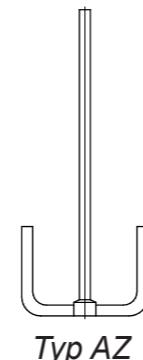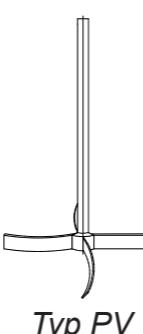
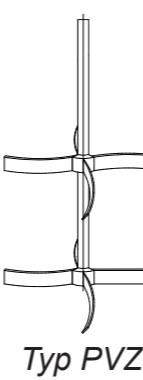
(3) aus V4A –

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(5) Deckel Typ D4303

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(1) Propeller stirrer with four blades and second stage

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- Following types are available: shaft diameter to flange size

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- or inner height of reactor and flange size

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- type, or article number, please have a look at page 82

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
(1) Anchor stirrer with two arms – (2) 8mm shaft – (3) out of V2A (PTFE coated) –

(4) for reactor type TRGN7168 – (5) lid type D4201 – (6) stirrer seal type RVFE 8/29

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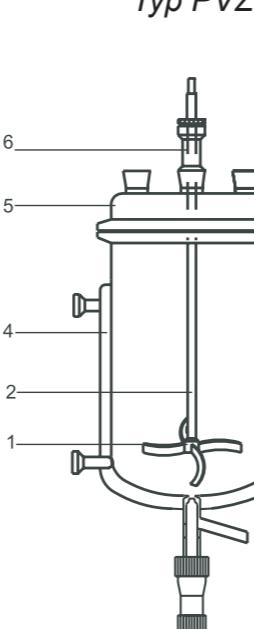
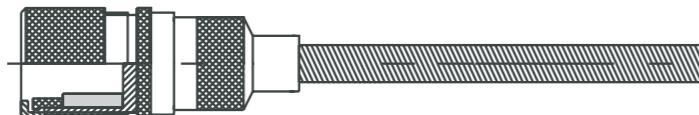
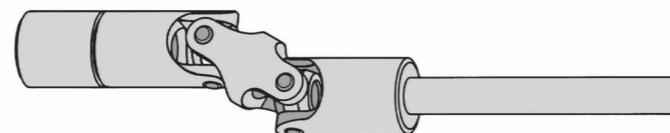
thermally insulated, with M16x1 threaded joints on both sides, matches to standard threads of circulation thermostats,
max. operating temperature 300° C

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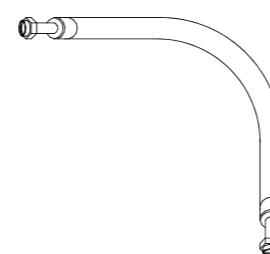
The KGW dynamic stirrer seal is suitable for both vacuum and excess pressure areas. It has two separate PTFE seals which can be adapted for excess pressure and vacuum applications by alternating of their mounting direction.

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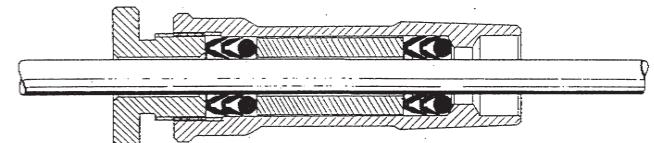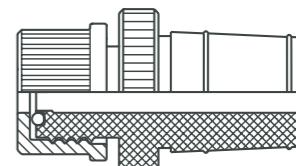
Vacuum and magnetic stirrer-seal coupling out of stainless steel 1.4401 CrNiMo, with permanent magnets, absolutely gas- and vacuum-tight acid- and alkali-resistant, maximum full speed 2.500 r.p.m.

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
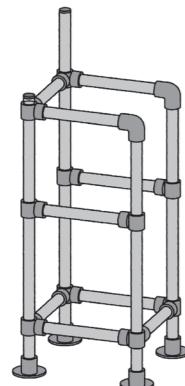
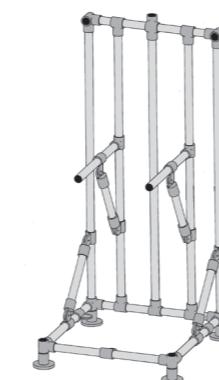
Advantages of the rack type Easy Frame

- fast changing of reactor
- only one person can easily install or remove the reactor
- no more loosening and dislocating of rack pipes while mounting the reactor, as necessary with ordinary frames
- reducing the risk of damage at flange or bottom outlet valve while mounting the reactor

The mobile rack type Easy Frame consists of:

- a stainless steel pipe frame with blue coated aluminum connectors
- a setting board out of aluminum and a insert with cushioning
- two adjustable flange holders for positioning and securing the reactor at its flange
- a stirrer machine rod and wheels

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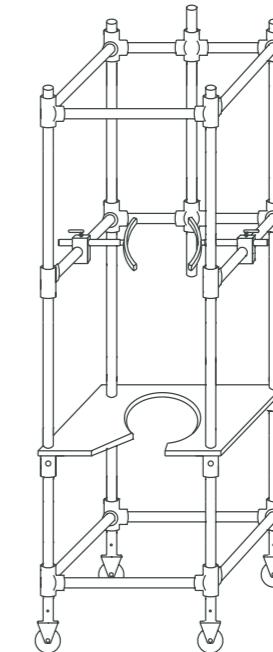The Reactor assembly for maintaining temperature consists of:
a heatable reactor with bottom outlet valve type „J“; metal adapters DN15 to male thread M16x1 for thermostat hoses and PTFE flat gaskets; a lid with three angular side sockets NS 29/32, one suitable centre socket and FEP coated O ring; a setting board; flange holders; V2A PTFE coated propeller stirrer; stirrer coupling RKG, dynamic stirrer seal RVFE and frame out of V2A tubes type Easy Frame.

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Further essential accessories: thermostat hoses, stirrer machine, thermostat.
On request we could also deliver an anchor stirrer.

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Reactor construction for maintaining temperature consists of the following:

a heatable reactor with 10 mm „J“ valve, metal adapter DN15 to M16x1 (PTFE sealants); a lid with three angular side sockets NS 29/32 and one suitable centre socket; an FEP coated O-ring; supporting device, V2A PTFE coated propeller stirrer, stirrer coupling RKG, dynamic stirrer seal RVFE and frame out of V2A tubes.

Further essential accessories: thermostat hoses, stirring machine, thermostat

On request we could also deliver an anchor stirrer.

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We would be happy to design your
individual assembly, you are welcome
to contact us.

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All columns produced by KGW-Isotherm are manufactured and tested to standard 12492 „Devices with vacuum insulation“. Each column is subjected to an individual test at a temperature of 200° C. All columns are manufactured with opposite viewing stripes as standard. The types listed in the tables are only a selection of our saleable columns.

All columns including the standard types are produced to the costumer's order. Therefore modifications to the ground glass joints are usually offered free of extra charge.

In addition to the types described here, other versions are available on request at any time. With the aid of CAD we are able to send you on request an offer with drawings.

Packed columns

with standard or spherical ground joints

Packed columns have excellent rectifying characteristics and ensure high throughputs for small operating volume and low pressure drops. A funnel for fixing of the packings is fused into the lower part of the column. This avoids bypassing at the edges of the column. When smaller packings are used, larger packings must first be inserted in the funnel to avoid that the smaller packings fall through the apertured funnel. The apertures in the funnel guarantees a free cross section of at least 50 percent of the base of ground joint.

Vigreux columns

with standard or spherical ground joints

Vigreux-columns are far more suitable for vacuum applications than packed columns, as the pressure difference is extremely low. Compared to packed columns, they have a very low rectifying effect.

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Technische Änderungen vorbehalten / All rights reserved for technical changes
All rights reserved for technical changes.